


hexyl 2-cyanoacrylate

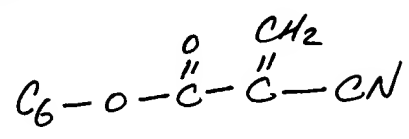

L3 1 S 3578-06-1
L4 21 S HEXYL AND CYANOACRYLATE

FILE 'CAPLUS, WPIDS, MEDLINE, EMBASE' ENTERED AT 16:27:11 ON 29 JUN 2002

FILE 'REGISTRY' ENTERED AT 16:27:58 ON 29 JUN 2002
SET SMARTSELECT ON

L5 SEL L3 1- CHEM : 6 TERMS
SET SMARTSELECT OFF

FILE 'CAPLUS, WPIDS, MEDLINE, EMBASE' ENTERED AT 16:27:59 ON 29 JUN 2002

L6 84 S L5/BI
L7 9304 S L6 OR CYANOACRYLATE# OR CYANACRYLATE#
L8 556 S L7 AND (AVM OR AVMS OR ARTERIOVENOUS MALFORMATION# OR OCCLUS
L9 274 S L7 AND (AVM OR AVMS OR ARTERIOVENOUS MALFORMATION#)
L10 6 S L9 AND PHOSPHORIC ACID
L11 4 DUP REM L10 (2 DUPLICATES REMOVED)
L12 7 S L8 AND PHOSPHORIC ACID
L13 1 S L12 NOT L10
L14 35 S (L8 AND (TANTALUM OR GOLD OR PLATINUM)) NOT (L10 OR L13)
L15 24 DUP REM L14 (11 DUPLICATES REMOVED)

=> d que 18

L3 1 SEA FILE=REGISTRY 3578-06-1
L5 SEL L3 1- CHEM : 6 TERMS
L6 84 SEA L5/BI
L7 9304 SEA L6 OR CYANOACRYLATE# OR CYANACRYLATE#
L8 556 SEA L7 AND (AVM OR AVMS OR ARTERIOVENOUS MALFORMATION# OR
OCCLUSION# OR OCCLUD?)

=> d 1-6 bib hit

L10 ANSWER 1 OF 6 CAPLUS COPYRIGHT 2002 ACS
AN 2002:123505 CAPLUS
DN 136:172869
TI Polymerizable **cyanoacrylate** compositions for therapeutic uses
IN Krall, Robert E.; Kerber, Charles W.; Knox, Kimberly
PA USA
SO U.S. Pat. Appl. Publ., 13 pp., Cont.-in-part of U.S. Ser. No. 577,115.
CODEN: USXXCO
DT Patent
LA English
FAN.CNT 2

	PATENT NO.	KIND	DATE	APPLICATION NO.	DATE
PI	US 2002018752	A1	20020214	US 2001-863825	20010523
PRAI	US 2000-577115	A2	20000523		

TI Polymerizable **cyanoacrylate** compositions for therapeutic uses
AB The present invention provides compns. that polymerizes upon contact with an anionic environment comprising (i) at least two polymerizable org. monomers and (ii) an oligomer of a polymerizable org. monomer, a plasticizer and an opacificant agent. Both polymerizable org. monomers are C1-18 alkyl **cyanoacrylates**. The compns. are useful for filling, occluding, partially filling or partially occluding an unfilled vol. or space in a mass in an anionic environment. The compn. are also useful for ablating diseased or undesired tissue, such as an **arteriovenous malformation** or a tumor, by cutting off the blood supply to the tissue. Also, the compns. are useful for controlled delivery of a therapeutic, chemotherapeutic or radiation delivery device to a desired location in the human body. For example, a formulation with monomer component was prepd. contg. n-hexyl **cyanoacrylate**, hydroquinone, p-methoxyphenol, and glacial acetic acid. The monomer component with a combination of Me **cyanoacrylate** and n-hexyl **cyanoacrylate** can be made. Twelve compounded **cyanoacrylates** were tested for conformal endovascular obliteration utility using a silicone aneurysm model: six based upon the 2-hexyl **cyanoacrylate**/methyl **cyanoacrylate** monomers, six based upon the 1-hexyl **cyanoacrylate**/methyl **cyanoacrylate** monomers. Additives consisted of various oils, gold for opacification, and polymn. retardants. All twelve compds. remained cohesive and conformed nicely to the outline of the aneurysm. Many of the mixts. based upon the 2-hexyl monomer exhibited delayed polymn., and could not be kept within the aneurysm lumen, even with adjacent balloon control of the infusion process. Four of the mixts. based upon the 1-hexyl monomer gave good cohesion, good conformation, remained within the aneurysm, and allowed some degree of angioplasty and remodeling of the arterial lumen by silicone balloon.
ST **cyanoacrylate** monomer oligomer polymn controlled drug delivery; blood vessel occlusion **cyanoacrylate** polymn
IT Polymerization
(anionic; polymerizable **cyanoacrylate** compns. for therapeutic uses)
IT Blood vessel, disease
(**arteriovenous malformation**; polymerizable **cyanoacrylate** compns. for endovascular occlusion in treatment of **arteriovenous malformation**)
IT Medical goods
(artificial venous valves; polymerizable **cyanoacrylate** compns. for tissue adhesion to medical surfaces)

IT Drug delivery systems
(controlled-release; polymerizable **cyanoacrylate** compns. for controlled drug delivery)

IT Fatty acids, biological studies
RL: THU (Therapeutic use); BIOL (Biological study); USES (Uses)
(esters, plasticizers; polymerizable **cyanoacrylate** compns. for therapeutic uses)

IT Sterility
(female, induction of; polymerizable **cyanoacrylate** compns. for administration to fallopian tubes in female sterilization)

IT Prosthetic materials and Prosthetics
(implants; polymerizable **cyanoacrylate** compns. for therapeutic uses)

IT Polymerization
(in situ; polymerizable **cyanoacrylate** compns. for therapeutic uses)

IT Castor oil
RL: THU (Therapeutic use); BIOL (Biological study); USES (Uses)
(iodinated; polymerizable **cyanoacrylate** compns. for therapeutic uses)

IT Uterus, neoplasm
(leiomyoma, inhibitors; polymerizable **cyanoacrylate** compns. for endovascular occlusion in tumor treatment)

IT Myoma
(leiomyoma, uterine, inhibitors; polymerizable **cyanoacrylate** compns. for endovascular occlusion in tumor treatment)

IT Blood vessel
(occlusion; polymerizable **cyanoacrylate** compns. for therapeutic uses)

IT Particle size
(of opacifiers; polymerizable **cyanoacrylate** compns. for therapeutic uses)

IT Magnetic materials
(particles, delivery of; polymerizable **cyanoacrylate** compns. for controlled drug delivery)

IT Drug delivery systems
(particles, magnetic; polymerizable **cyanoacrylate** compns. for controlled drug delivery)

IT Oviduct
(polymerizable **cyanoacrylate** compns. for administration to fallopian tubes in female sterilization)

IT Chemotherapy
Gene therapy
Radiotherapy
(polymerizable **cyanoacrylate** compns. for controlled drug delivery)

IT Antitumor agents
(polymerizable **cyanoacrylate** compns. for endovascular occlusion in tumor treatment)

IT Aneurysm
Human
Opacifiers
Plasticizers
Polymerization inhibitors
(polymerizable **cyanoacrylate** compns. for therapeutic uses)

IT Monomers
Oligomers
RL: RCT (Reactant); THU (Therapeutic use); BIOL (Biological study); RACT (Reactant or reagent); USES (Uses)

(polymerizable **cyanoacrylate** compns. for therapeutic uses)

IT Acids, biological studies
 RL: THU (Therapeutic use); BIOL (Biological study); USES (Uses)
 (polymn. inhibitors; polymerizable **cyanoacrylate** compns. for therapeutic uses)

IT Polymerization
 (radical; polymerizable **cyanoacrylate** compns. for therapeutic uses)

IT Medical goods
 (stents; polymerizable **cyanoacrylate** compns. for tissue adhesion to medical surfaces)

IT Antitumor agents
 (uterus leiomyoma; polymerizable **cyanoacrylate** compns. for endovascular occlusion in tumor treatment)

IT Heart
 (valve, artificial; polymerizable **cyanoacrylate** compns. for tissue adhesion to medical surfaces)

IT Vein
 (valves, artificial; polymerizable **cyanoacrylate** compns. for tissue adhesion to medical surfaces)

IT 7440-05-3, Palladium, biological studies 7440-06-4, Platinum, biological studies 7440-25-7, Tantalum, biological studies 7440-32-6, Titanium, biological studies 7440-57-5, Gold, biological studies
 RL: THU (Therapeutic use); BIOL (Biological study); USES (Uses)
 (opacifier; polymerizable **cyanoacrylate** compns. for therapeutic uses)

IT 57-10-3D, Palmitic acid, esters 57-11-4D, Stearic acid, esters 112-80-1D, Oleic acid, esters 124-06-1, Ethyl myristate 143-07-7D, Lauric acid, esters 544-63-8D, Myristic acid, esters
 RL: THU (Therapeutic use); BIOL (Biological study); USES (Uses)
 (plasticizer; polymerizable **cyanoacrylate** compns. for therapeutic uses)

IT 137-05-3, Methyl **cyanoacrylate** 1069-55-2, Isobutyl **cyanoacrylate** 3578-06-1 6606-65-1 6701-17-3, n-Octyl 2-**cyanoacrylate** 15721-32-1, 2-Ethylhexyl 2-**cyanoacrylate** 15802-18-3D, Cyanoacrylic acid, alkyl esters 398147-86-9
 RL: RCT (Reactant); THU (Therapeutic use); BIOL (Biological study); RACT (Reactant or reagent); USES (Uses)
 (polymerizable **cyanoacrylate** compns. for therapeutic uses)

IT 123-31-9, Hydroquinone, biological studies 150-76-5, p-Methoxyphenol
 RL: THU (Therapeutic use); BIOL (Biological study); USES (Uses)
 (polymerizable **cyanoacrylate** compns. for therapeutic uses)

IT 64-19-7, Acetic acid, biological studies 7664-38-2, **Phosphoric acid**, biological studies
 RL: THU (Therapeutic use); BIOL (Biological study); USES (Uses)
 (polymn. inhibitor; polymerizable **cyanoacrylate** compns. for therapeutic uses)

L10 ANSWER 2 OF 6 CAPLUS COPYRIGHT 2002 ACS

AN 2001:868200 CAPLUS

DN 136:11252

TI Polymerizable compositions as filling materials and methods of their use

IN Krall, Robert E.; Kerber, Charles W.; Knox, Kimberly

PA Provasis Therapeutics, Inc., USA

SO PCT Int. Appl., 47 pp.

CODEN: PIXXD2

DT Patent

LA English

FAN.CNT 2

	PATENT NO.	KIND	DATE	APPLICATION NO.	DATE
PI	WO 2001089501	A1	20011129	WO 2001-US16638	20010523
	W:	AE, AG, AL, AM, AT, AU, AZ, BA, BB, BG, BR, BY, BZ, CA, CH, CN, CO, CR, CU, CZ, DE, DK, DM, DZ, EC, EE, ES, FI, GB, GD, GE, GH, GM, HR, HU, ID, IL, IN, IS, JP, KE, KG, KP, KR, KZ, LC, LK, LR, LS, LT, LU, LV, MA, MD, MG, MK, MN, MW, MX, MZ, NO, NZ, PL, PT, RO, RU, SD, SE, SG, SI, SK, SL, TJ, TM, TR, TT, TZ, UA, UG, UZ, VN, YU, ZA, ZW, AM, AZ, BY, KG, KZ, MD, RU, TJ, TM			
	RW:	GH, GM, KE, LS, MW, MZ, SD, SL, SZ, TZ, UG, ZW, AT, BE, CH, CY, DE, DK, ES, FI, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE, TR, BF, BJ, CF, CG, CI, CM, GA, GN, GW, ML, MR, NE, SN, TD, TG			

PRAI US 2000-577115 A 20000523

RE.CNT 10 THERE ARE 10 CITED REFERENCES AVAILABLE FOR THIS RECORD

ALL CITATIONS AVAILABLE IN THE RE FORMAT

AB The present invention provides compns. comprising a first component and a second component, wherein the first component includes at least two polymerizable org. monomers, and wherein the second component includes an oligomer of a polymerizable org. monomer, a plasticizer and an opacificant agent, wherein said compn. polymerizes upon contact with an anionic environment. The compns. are useful for filling, occluding, partially filling or partially occluding an unfilled vol. or space in a mass in an anionic environment. The compn. are also useful for ablating diseased or undesired tissue by cutting off the blood supply to the tissue. A first component contained 2-hexyl **cyanoacrylate** 1250, hydroquinone 0.0764, p-methoxyphenol 0.0874, and **phosphoric acid** 0.1693 g. A second component was prepd. by mixing 2.0 g of oligo(2-hexyl **cyanoacrylate**) with 100 g of powd. gold, then 1.020 g of this blended material was mixed with 500 mg of Et myristate. Comparison of catheter adhesion force demonstrated that 2-hexyl **cyanoacrylate** compn. had significantly lower adhesion to the catheter than the controls contg. Bu **cyanoacrylate**.

ST polymerizable filling material alkyl **cyanoacrylate**

IT Blood vessel, disease

(arteriovenous malformation; polymerizable compns. as filling materials and methods of their use)

IT 15802-18-3DP, Cyanoacrylic acid, alkyl derivs. 25154-80-7P, Poly(butyl **cyanoacrylate**) 26809-38-1P, Poly(iso-butyl **cyanoacrylate**) 26877-34-9P, Poly(octyl **cyanoacrylate**) 26877-39-4P, Poly(2-Hexyl **cyanoacrylate**)

RL: DEV (Device component use); SPN (Synthetic preparation); THU (Therapeutic use); BIOL (Biological study); PREP (Preparation); USES (Uses)

(polymerizable compns. as filling materials and methods of their use)

IT 57-10-3D, Palmitic acid, esters 57-11-4D, Stearic acid, esters 64-19-7, Acetic acid, uses 112-80-1D, Oleic acid, esters 124-06-1, Ethyl myristate 143-07-7D, Lauric acid, esters 544-63-8D, Myristic acid, esters 7440-05-3, Palladium, uses 7440-06-4, Platinum, uses 7440-25-7, Tantalum, uses 7440-32-6, Titanium, uses 7440-57-5, Gold, uses 7664-38-2, **Phosphoric acid**, uses

RL: NUU (Other use, unclassified); POF (Polymer in formulation); USES (Uses)

(polymerizable compns. as filling materials and methods of their use)

L10 ANSWER 3 OF 6 CAPLUS COPYRIGHT 2002 ACS

AN 2000:534949 CAPLUS

DN 133:140311

TI **Cyanoacrylates** comprising inhibitors and an opacifying agent as

adhesives

IN Krall, Robert E.; Kerber, Charles W.; Knox, Kimberly

PA Prohold Medical Technologies, Inc., USA

SO PCT Int. Appl., 70 pp.

CODEN: PIXXD2

DT Patent

LA English

FAN.CNT 1

	PATENT NO.	KIND	DATE	APPLICATION NO.	DATE
PI	WO 2000044287	A1	20000803	WO 2000-US2262	20000128
	W:	AE, AL, AM, AT, AU, AZ, BA, BB, BG, BR, BY, CA, CH, CN, CR, CU, CZ, DE, DK, DM, EE, ES, FI, GB, GD, GE, GH, GM, HR, HU, ID, IL, IN, IS, JP, KE, KG, KP, KR, KZ, LC, LK, LR, LS, LT, LU, LV, MA, MD, MG, MK, MN, MW, MX, NO, NZ, PL, PT, RO, RU, SD, SE, SG, SI, SK, SL, TJ, TM, TR, TT, TZ, UA, UG, US, UZ, VN, YU, ZA, ZW, AM, AZ, BY, KG, KZ, MD, RU, TJ, TM			
	RW:	GH, GM, KE, LS, MW, SD, SL, SZ, TZ, UG, ZW, AT, BE, CH, CY, DE, DK, ES, FI, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE, BF, BJ, CF, CG, CI, CM, GA, GN, GW, ML, MR, NE, SN, TD, TG			
	EP 1154723	A1	20011121	EP 2000-904626	20000128
	R:	AT, BE, CH, DE, DK, ES, FR, GB, GR, IT, LI, LU, NL, SE, MC, PT, IE, SI, LT, LV, FI, RO			
PRAI	US 1999-241368	A2	19990129		
	WO 2000-US2262	W	20000128		

RE.CNT 6 THERE ARE 6 CITED REFERENCES AVAILABLE FOR THIS RECORD

ALL CITATIONS AVAILABLE IN THE RE FORMAT

TI **Cyanoacrylates** comprising inhibitors and an opacifying agent as adhesives

AB A compn. comprising of a monomer component comprised of an alkyl **cyanoacrylate** and at least one inhibitor, and a second component comprised of a resultant aggregate structure formed from an alkyl **cyanoacrylate** monomer, an alkyl esterified fatty acid and an opacificant agent where said compn. forms a resultant aggregate structure when said compn. contacts an anionic environment. The compn. is useful for filling an existing space, e.g., the lumen of a blood vessel, a space created by a transiently placed external device, e.g., a catheter or like device, a space created by a procedure, e.g., an excision or implantation of an object, e.g., a stent. The compn. is also useful for adhering tissue to tissue, or adhering tissue to a device. The compn. has the property of polymg. when it comes in contact with an anionic environment, or when it is deployed in situ in an existing space. Thus, 2-hexyl **cyanoacrylate** (I) was prepd. by the reaction of paraformaldehyde with 2-hexyl cyanoacetate. An adhesive formulation contained I 6.8964, hydroquinone 0.000694, p-methoxyphenol 0.000704, **phosphoric acid** 0.001726 mol. The formulation was an effective embolic agent in the treatment of a patient after acute hemorrhage of a right parieto-occipital **arteriovenous malformation**.

ST **cyanoacrylate** inhibitor opacifier tissue adhesive

IT Embolism

(agents for; **cyanoacrylates** comprising inhibitors and opacifying agent as adhesives)

IT Adhesives

Adhesives

(biol. tissue; **cyanoacrylates** comprising inhibitors and opacifying agent as adhesives)

IT Opacifiers

(**cyanoacrylates** comprising inhibitors and opacifying agent as adhesives)

IT Fatty acids, biological studies
 RL: POF (Polymer in formulation); THU (Therapeutic use); BIOL (Biological study); USES (Uses)
 (cyanoacrylates comprising inhibitors and opacifying agent as adhesives)

IT Liquids
 (oils, halogenated; cyanoacrylates comprising inhibitors and opacifying agent as adhesives)

IT Medical goods
 Medical goods
 (tissue adhesives; cyanoacrylates comprising inhibitors and opacifying agent as adhesives)

IT 3578-06-1P, 2-Hexyl cyanoacrylate
 RL: BAC (Biological activity or effector, except adverse); BSU (Biological study, unclassified); POF (Polymer in formulation); SPN (Synthetic preparation); THU (Therapeutic use); BIOL (Biological study); PREP (Preparation); USES (Uses)
 (cyanoacrylates comprising inhibitors and opacifying agent as adhesives)

IT 57-10-3D, Palmitic acid, alkyl derivs. 57-11-4D, Stearic acid, alkyl derivs. 123-31-9, 1,4-Benzenediol, biological studies 124-06-1, Ethyl myristate 143-07-7D, Lauric acid, alkyl derivs. 150-76-5, p-Methoxyphenol 7440-06-4, Platinum, biological studies 7440-25-7, Tantalum, biological studies 7440-32-6, Titanium, biological studies 7440-33-7, Tungsten, biological studies 7440-57-5, Gold, biological studies 7664-38-2, Phosphoric acid, biological studies 7727-43-7, Barium sulfate 8008-53-5, Ethiodol
 RL: POF (Polymer in formulation); THU (Therapeutic use); BIOL (Biological study); USES (Uses)
 (cyanoacrylates comprising inhibitors and opacifying agent as adhesives)

IT 372-09-8, Cyanoacetic acid 626-93-7, 2-Hexanol 30525-89-4, Paraformaldehyde
 RL: RCT (Reactant); RACT (Reactant or reagent)
 (cyanoacrylates comprising inhibitors and opacifying agent as adhesives)

IT 286931-66-6P
 RL: RCT (Reactant); SPN (Synthetic preparation); PREP (Preparation); RACT (Reactant or reagent)
 (cyanoacrylates comprising inhibitors and opacifying agent as adhesives)

L10 ANSWER 4 OF 6 CAPLUS COPYRIGHT 2002 ACS

AN 2000:172844 CAPLUS

DN 132:212712

TI Composition comprising 2-hexyl cyanoacrylate and gold for creating vascular occlusions

IN Krall, Robert E.; Kerber, Charles W.; Knox, Kimberly

PA Prohold Medical Technologies, Inc., USA

SO U.S., 3 pp.

CODEN: USXXAM

DT Patent

LA English

FAN.CNT 1

	PATENT NO.	KIND	DATE	APPLICATION NO.	DATE
	-----	---	-----	-----	-----
PI	US 6037366	A	20000314	US 1998-151621	19980911
PRAI	US 1997-58510P	P	19970911		
RE.CNT	6				

THERE ARE 6 CITED REFERENCES AVAILABLE FOR THIS RECORD

ALL CITATIONS AVAILABLE IN THE RE FORMAT

TI Composition comprising 2-hexyl **cyanoacrylate** and gold for creating vascular occlusions

AB A compn. including 2-hexyl **cyanoacrylate** and gold is useful in treating **arteriovenous malformations (AVMs)** and other body lumens to be blocked. A compn. comprised 2-hexyl **cyanoacrylate** 999.550, hydroquinone 100, methoxyphenol 100, **phosphoric acid** 250 ppm in the first part; and pure gold 1.0000, pure Et myristate 0.5000, and FMS (a specially prepd. polymer of 2-hexyl **cyanoacrylate**) 0.0200 g in the second part.

ST vascular occlusion **cyanoacrylate** gold **arteriovenous malformation**

IT Blood vessel, disease
(occlusion; compn. comprising 2-hexyl **cyanoacrylate** and gold for creating vascular occlusions)

IT 7440-57-5, Gold, biological studies
RL: BAC (Biological activity or effector, except adverse); BSU (Biological study, unclassified); THU (Therapeutic use); BIOL (Biological study); USES (Uses)
(compn. comprising 2-hexyl **cyanoacrylate** and gold for creating vascular occlusions)

IT 123-31-9, Hydroquinone, biological studies 124-06-1, Ethyl myristate **3578-06-1**, 2-Hexyl **cyanoacrylate** 7664-38-2, **Phosphoric acid**, biological studies 26638-03-9, Methoxyphenol
RL: THU (Therapeutic use); BIOL (Biological study); USES (Uses)
(compn. comprising 2-hexyl **cyanoacrylate** and gold for creating vascular occlusions)

L10 ANSWER 5 OF 6 WPIDS (C) 2002 THOMSON DERWENT

AN 2002-240636 [29] WPIDS

DNC C2002-072360

TI Polymerizable composition for filling or occluding unfilled volume or space in human or animal body, contains first component comprising polymerizable organic monomers, and second component containing organic oligomer.

DC A14 A96 B04 D22

IN KERBER, C W; KNOX, K; KRALL, R E

PA (PROV-N) PROVAVIS THERAPEUTICS INC

CYC 95

PI WO 2001089501 A1 20011129 (200229)* EN 47p
RW: AT BE CH CY DE DK EA ES FI FR GB GH GM GR IE IT KE LS LU MC MW MZ
NL OA PT SD SE SL SZ TR TZ UG ZW
W: AE AG AL AM AT AU AZ BA BB BG BR BY BZ CA CH CN CO CR CU CZ DE DK
DM DZ EC EE ES FI GB GD GE GH GM HR HU ID IL IN IS JP KE KG KP KR
KZ LC LK LR LS LT LU LV MA MD MG MK MN MW MX MZ NO NZ PL PT RO RU
SD SE SG SI SK SL TJ TM TR TT TZ UA UG UZ VN YU ZA ZW
AU 2001064835 A 20011203 (200229)

ADT WO 2001089501 A1 WO 2001-US16638 20010523; AU 2001064835 A AU 2001-64835 20010523

FDT AU 2001064835 A Based on WO 200189501

PRAI US 2000-577115 20000523

AB WO 200189501 A UPAB: 20020508
NOVELTY - A polymerizable composition consists of first and second components. The first component comprises at least 2 polymerizable organic monomers, and the second component contains an oligomer of polymerizable organic monomer, a plasticizer, and an opacificant agent. The composition polymerizes upon contact with an anionic environment.
USE - The polymerizable composition is used for filling or occluding

an unfilled volume or space e.g. the lumen of blood vessel, the sac of an aneurysm, a space created by transiently placed external device such as catheter, or a space created by a procedure such as excision or implantation of stent. It is also useful for ablating diseased or undesired tissue e.g. **arteriovenous malformation**, and for treating tumor and uterine leiomyoma by blocking or cutting off the blood supply to the tissue or organ. The composition is also used for adhering tissue to tissue, or tissue to a device (claimed). It is also useful as an embolic agent that selectively creates an embolic blockage in blood vessel lumen, duct, fistula, or other body passageways. The composition can be used for aortopulmonary closure; treatment of artery pseudoaneurysm; hepatic artery vascular occlusion and temporary vascular occlusion during co-administration of cytotoxic drugs; and for creating tubal occlusion, fallopian tube occlusion, vas deferens occlusion, and urinary occlusion.

ADVANTAGE - The polymerizable composition is radiopaque, and the heat released during its polymerization does not adversely affect heat-sensitive surrounding tissues e.g. brain tissues. The composition and its biodegradation products are non-histotoxic and non-cytotoxic.

Dwg.0/0

TECH

UPTX: 20020508

TECHNOLOGY FOCUS - POLYMERS - Preferred Components: The first component comprises 1-500 (200-300) ppm polymerization inhibitor(s) (preferably acetic or **phosphoric acid**) to inhibit free radical or anionic polymerization. The second component may contain halogenated oil (preferably iodinated castor oil).

Preferred Compounds: The polymerizable organic monomers are 1-18C alkyl **cyanoacrylates**, preferably methyl **cyanoacrylates**, n-butyl **cyanoacrylates**, isobutyl **cyanoacrylates**, n-hexyl **cyanoacrylates**, 2-hexyl **cyanoacrylates**, n-octyl **cyanoacrylates**, or 2-ethylhexyl **cyanoacrylates**. The plasticizer is an esterified fatty acid e.g. laurate, palmitate, oleate, myristate, or stearate, preferably ethyl myristate.

TECHNOLOGY FOCUS - INORGANIC CHEMISTRY - Preferred Compound: The opacificant agent is a metal e.g. gold, platinum, palladium, tantalum, and/or titanium, or its alloy, or preferably gold in fine powder form having particle diameter of not more than 7 (preferably not more than 1) microns.

L10 ANSWER 6 OF 6 WPIDS (C) 2002 THOMSON DERWENT

AN 2000-246196 [21] WPIDS

DNC C2000-074484

TI Composition comprising a 2-hexyl **cyanoacrylate** monomer and gold in a polymer of 2-hexylcyanoacrylate, useful for placing in a body lumen to create vascular occlusion.

DC A14 A96 B05

IN KERBER, C W; KNOX, K; KRALL, R E

PA (PROH-N) PROHOLD MEDICAL TECHNOLOGIES INC

CYC 1

PI US 6037366 A 20000314 (200021)* 3p

ADT US 6037366 A Provisional US 1997-58510P 19970911, US 1998-151621 19980911

PRAI US 1997-58510P 19970911; US 1998-151621 19980911

TI Composition comprising a 2-hexyl **cyanoacrylate** monomer and gold in a polymer of 2-hexylcyanoacrylate, useful for placing in a body lumen to create vascular occlusion.

AB US 6037366 A UPAB: 20000502

NOVELTY - A composition for placing in a body lumen to create vascular occlusion comprises a 2-hexyl **cyanoacrylate** monomer and gold in

a polymer of 2-hexylcyanoacrylate.

DETAILED DESCRIPTION - A composition for creating vascular occlusions comprises a mixture of:

(a) 2-hexyl **cyanoacrylate**, hydroquinone, p-methoxyphenol and **phosphoric acid**; and

(b) gold metal powder, ethyl myristate, and a sterilized polymer of 2-hexylcyanoacrylate in weak aqueous bicarbonate solution.

ACTIVITY - Cytostatic; Antiarteriosclerotic; Vasotropic; Cerebroprotective.

MECHANISM OF ACTION - None given.

No biological data is given.

USE - The composition is useful for creating vascular occlusions and for treating **arteriovenous malformations** and tumors (particularly neurological).

ADVANTAGE - The composition is especially useful to treat vascular tumors in the brain and brain stem, both of which are difficult to access and which are susceptible to cytotoxicity and heat.
Dwg.0/0

TECH UPTX: 20000502

TECHNOLOGY FOCUS - PHARMACEUTICALS - Preferred Composition: Component (a) comprises 100 parts per million (ppm) hydroquinone, 100 ppm p-methoxyphenol, 250 ppm **phosphoric acid** and the remainder 2-hexyl **cyanoacrylate**. Component (b) comprises 65 wt.% gold, 30 wt.% ethyl myristate and the remainder sterilized polymer of 2-hexylcyanoacrylate in weak aqueous bicarbonate solution, optionally with sulfur dioxide as a stabilizer.

TT TT: COMPOSITION COMPRISE HEXYL **CYANOACRYLATE** MONOMER GOLD
POLYMER USEFUL PLACE BODY LUMEN VASCULAR OCCLUDE.

=> d 1-4 bib ab

L11 ANSWER 1 OF 4 CAPLUS COPYRIGHT 2002 ACS
AN 2002:123505 CAPLUS
DN 136:172869
TI Polymerizable **cyanoacrylate** compositions for therapeutic uses
IN Krall, Robert E.; Kerber, Charles W.; Knox, Kimberly
PA USA
SO U.S. Pat. Appl. Publ., 13 pp., Cont.-in-part of U.S. Ser. No. 577,115.
CODEN: USXXCO
DT Patent
LA English
FAN.CNT 2

	PATENT NO.	KIND	DATE	APPLICATION NO.	DATE
PI	US 2002018752	A1	20020214	US 2001-863825	20010523
PRAI	US 2000-577115	A2	20000523		

AB The present invention provides compns. that polymerizes upon contact with an anionic environment comprising (i) at least two polymerizable org. monomers and (ii) an oligomer of a polymerizable org. monomer, a plasticizer and an opacificant agent. Both polymerizable org. monomers are C1-18 alkyl **cyanoacrylates**. The compns. are useful for filling, occluding, partially filling or partially occluding an unfilled vol. or space in a mass in an anionic environment. The compn. are also useful for ablating diseased or undesired tissue, such as an **arteriovenous malformation** or a tumor, by cutting off the blood supply to the tissue. Also, the compns. are useful for controlled delivery of a therapeutic, chemotherapeutic or radiation delivery device to a desired location in the human body. For example, a formulation with monomer component was prepd. contg. n-hexyl **cyanoacrylate**, hydroquinone, p-methoxyphenol, and glacial acetic acid. The monomer component with a combination of Me **cyanoacrylate** and n-hexyl **cyanoacrylate** can be made. Twelve compounded **cyanoacrylates** were tested for conformal endovascular obliteration utility using a silicone aneurysm model: six based upon the 2-hexyl **cyanoacrylate**/methyl **cyanoacrylate** monomers, six based upon the 1-hexyl **cyanoacrylate**/methyl **cyanoacrylate** monomers. Additives consisted of various oils, gold for opacification, and polymn. retardants. All twelve compds. remained cohesive and conformed nicely to the outline of the aneurysm. Many of the mixts. based upon the 2-hexyl monomer exhibited delayed polymn., and could not be kept within the aneurysm lumen, even with adjacent balloon control of the infusion process. Four of the mixts. based upon the 1-hexyl monomer gave good cohesion, good conformation, remained within the aneurysm, and allowed some degree of angioplasty and remodeling of the arterial lumen by silicone balloon.

L11 ANSWER 2 OF 4 CAPLUS COPYRIGHT 2002 ACS DUPLICATE 1
AN 2001:868200 CAPLUS
DN 136:11252
TI Polymerizable compositions as filling materials and methods of their use
IN Krall, Robert E.; Kerber, Charles W.; Knox, Kimberly
PA Provasis Therapeutics, Inc., USA
SO PCT Int. Appl., 47 pp.
CODEN: PIXXD2
DT Patent
LA English
FAN.CNT 2

	PATENT NO.	KIND	DATE	APPLICATION NO.	DATE
--	------------	------	------	-----------------	------

 PI WO 2001089501 A1 20011129 WO 2001-US16638 20010523
 W: AE, AG, AL, AM, AT, AU, AZ, BA, BB, BG, BR, BY, BZ, CA, CH, CN,
 CO, CR, CU, CZ, DE, DK, DM, DZ, EC, EE, ES, FI, GB, GD, GE, GH,
 GM, HR, HU, ID, IL, IN, IS, JP, KE, KG, KP, KR, KZ, LC, LK, LR,
 LS, LT, LU, LV, MA, MD, MG, MK, MN, MW, MX, MZ, NO, NZ, PL, PT,
 RO, RU, SD, SE, SG, SI, SK, SL, TJ, TM, TR, TT, TZ, UA, UG, UZ,
 VN, YU, ZA, ZW, AM, AZ, BY, KG, KZ, MD, RU, TJ, TM
 RW: GH, GM, KE, LS, MW, MZ, SD, SL, SZ, TZ, UG, ZW, AT, BE, CH, CY,
 DE, DK, ES, FI, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE, TR, BF,
 BJ, CF, CG, CI, CM, GA, GN, GW, ML, MR, NE, SN, TD, TG

PRAI US 2000-577115 A 20000523

AB The present invention provides compns. comprising a first component and a second component, wherein the first component includes at least two polymerizable org. monomers, and wherein the second component includes an oligomer of a polymerizable org. monomer, a plasticizer and an opacificant agent, wherein said compn. polymerizes upon contact with an anionic environment. The compns. are useful for filling, occluding, partially filling or partially occluding an unfilled vol. or space in a mass in an anionic environment. The compn. are also useful for ablating diseased or undesired tissue by cutting off the blood supply to the tissue. A first component contained 2-hexyl **cyanoacrylate** 1250, hydroquinone 0.0764, p-methoxyphenol 0.0874, and **phosphoric acid** 0.1693 g. A second component was prepd. by mixing 2.0 g of oligo(2-hexyl **cyanoacrylate**) with 100 g of powd. gold, then 1.020 g of this blended material was mixed with 500 mg of Et myristate. Comparison of catheter adhesion force demonstrated that 2-hexyl **cyanoacrylate** compn. had significantly lower adhesion to the catheter than the controls contg. Bu **cyanoacrylate**.

RE.CNT 10 THERE ARE 10 CITED REFERENCES AVAILABLE FOR THIS RECORD
 ALL CITATIONS AVAILABLE IN THE RE FORMAT

L11 ANSWER 3 OF 4 CAPLUS COPYRIGHT 2002 ACS DUPLICATE 2

AN 2000:172844 CAPLUS

DN 132:212712

TI Composition comprising 2-hexyl **cyanoacrylate** and gold for creating vascular occlusions

IN Krall, Robert E.; Kerber, Charles W.; Knox, Kimberly

PA Prohold Medical Technologies, Inc., USA

SO U.S., 3 pp.

CODEN: USXXAM

DT Patent

LA English

FAN.CNT 1

	PATENT NO.	KIND	DATE	APPLICATION NO.	DATE
PI	US 6037366	A	20000314	US 1998-151621	19980911
PRAI	US 1997-58510P	P	19970911		

AB A compn. including 2-hexyl **cyanoacrylate** and gold is useful in treating **arteriovenous malformations (AVMs)** and other body lumens to be blocked. A compn. comprised 2-hexyl **cyanoacrylate** 999.550, hydroquinone 100, methoxyphenol 100, **phosphoric acid** 250 ppm in the first part; and pure gold 1.0000, pure Et myristate 0.5000, and FMS (a specially prepd. polymer of 2-hexyl **cyanoacrylate**) 0.0200 g in the second part.

RE.CNT 6 THERE ARE 6 CITED REFERENCES AVAILABLE FOR THIS RECORD
 ALL CITATIONS AVAILABLE IN THE RE FORMAT

L11 ANSWER 4 OF 4 CAPLUS COPYRIGHT 2002 ACS

AN 2000:534949 CAPLUS
 DN 133:140311
 TI **Cyanoacrylates** comprising inhibitors and an opacifying agent as adhesives
 IN Krall, Robert E.; Kerber, Charles W.; Knox, Kimberly
 PA Prohold Medical Technologies, Inc., USA
 SO PCT Int. Appl., 70 pp.
 CODEN: PIXXD2
 DT Patent
 LA English
 FAN.CNT 1

	PATENT NO.	KIND	DATE	APPLICATION NO.	DATE
PI	WO 2000044287	A1	20000803	WO 2000-US2262	20000128
	W:	AE, AL, AM, AT, AU, AZ, BA, BB, BG, BR, BY, CA, CH, CN, CR, CU, CZ, DE, DK, DM, EE, ES, FI, GB, GD, GE, GH, GM, HR, HU, ID, IL, IN, IS, JP, KE, KG, KP, KR, KZ, LC, LK, LR, LS, LT, LU, LV, MA, MD, MG, MK, MN, MW, MX, NO, NZ, PL, PT, RO, RU, SD, SE, SG, SI, SK, SL, TJ, TM, TR, TT, TZ, UA, UG, US, UZ, VN, YU, ZA, ZW, AM, AZ, BY, KG, KZ, MD, RU, TJ, TM			
	RW:	GH, GM, KE, LS, MW, SD, SL, SZ, TZ, UG, ZW, AT, BE, CH, CY, DE, DK, ES, FI, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE, BF, BJ, CF, CG, CI, CM, GA, GN, GW, ML, MR, NE, SN, TD, TG			
	EP 1154723	A1	20011121	EP 2000-904626	20000128
	R:	AT, BE, CH, DE, DK, ES, FR, GB, GR, IT, LI, LU, NL, SE, MC, PT, IE, SI, LT, LV, FI, RO			
PRAI	US 1999-241368	A2	19990129		
	WO 2000-US2262	W	20000128		
AB	<p>A compn. comprising of a monomer component comprised of an alkyl cyanoacrylate and at least one inhibitor, and a second component comprised of a resultant aggregate structure formed from an alkyl cyanoacrylate monomer, an alkyl esterified fatty acid and an opacificant agent where said compn. forms a resultant aggregate structure when said compn. contacts an anionic environment. The compn. is useful for filling an existing space, e.g., the lumen of a blood vessel, a space created by a transiently placed external device, e.g., a catheter or like device, a space created by a procedure, e.g., an excision or implantation of an object, e.g., a stent. The compn. is also useful for adhering tissue to tissue, or adhering tissue to a device. The compn. has the property of polymg. when it comes in contact with an anionic environment, or when it is deployed in situ in an existing space. Thus, 2-hexyl cyanoacrylate (I) was prepd. by the reaction of paraformaldehyde with 2-hexyl cyanoacetate. An adhesive formulation contained I 6.8964, hydroquinone 0.000694, p-methoxyphenol 0.000704, phosphoric acid 0.001726 mol. The formulation was an effective embolic agent in the treatment of a patient after acute hemorrhage of a right parieto-occipital arteriovenous malformation.</p>				

RE.CNT 6 THERE ARE 6 CITED REFERENCES AVAILABLE FOR THIS RECORD
 ALL CITATIONS AVAILABLE IN THE RE FORMAT

L13 ANSWER 1 OF 1 WPIDS (C) 2002 THOMSON DERWENT

AN 2000-543304 [49] WPIDS

DNN N2000-401955 DNC C2000-161579

TI Composition which includes an alkyl **cyanoacrylate** monomer is used e.g. for adhering skin tissue.

DC A14 A96 D22 E14 E36 G03 P31 P34

IN KERBER, C W; KNOX, K; KRALL, R E

PA (PROH-N) PROHOLD MEDICAL TECHNOLOGIES INC

CYC 91

PI WO 2000044287 A1 20000803 (200049)* EN 65p

RW: AT BE CH CY DE DK EA ES FI FR GB GH GM GR IE IT KE LS LU MC MW NL
OA PT SD SE SL SZ TZ UG ZW

W: AE AL AM AT AU AZ BA BB BG BR BY CA CH CN CR CU CZ DE DK DM EE ES
FI GB GD GE GH GM HR HU ID IL IN IS JP KE KG KP KR KZ LC LK LR LS
LT LU LV MA MD MG MK MN MW MX NO NZ PL PT RO RU SD SE SG SI SK SL
TJ TM TR TT TZ UA UG US UZ VN YU ZA ZW

AU 2000026354 A 20000818 (200057)

EP 1154723 A1 20011121 (200176) EN

R: AL AT BE CH CY DE DK ES FI FR GB GR IE IT LI LT LU LV MC MK NL PT
RO SE SI

ADT WO 2000044287 A1 WO 2000-US2262 20000128; AU 2000026354 A AU 2000-26354
20000128; EP 1154723 A1 EP 2000-904626 20000128, WO 2000-US2262 20000128

FDT AU 2000026354 A Based on WO 200044287; EP 1154723 A1 Based on WO 200044287

PRAI US 1999-241368 19990129

TI Composition which includes an alkyl **cyanoacrylate** monomer is used e.g. for adhering skin tissue.

AB WO 200044287 A UPAB: 20001006

NOVELTY - A composition which forms an aggregate structure when contacted with an anionic environment comprises an alkyl **cyanoacrylate** monomer with at least one inhibitor, as well as a component formed from an alkyl **cyanoacrylate** monomer, an alkyl esterified fatty acid and an opacifying agent.

DETAILED DESCRIPTION - INDEPENDENT CLAIM is included for a method for filling, **occluding**, partially filling or partially **occluding** an unfilled volume or space in a mass in an anionic environment by using the claimed composition.

USE - The composition is useful for filling an existing space, such as the lumen of a blood vessel, a space created by a transiently placed external device, such as a catheter, or a space created by an excision or the implantation of an object, such as a stent. The composition is also useful for adhering tissue to tissue or tissue to a device.

ADVANTAGE - The composition has desirable cohesion, stability, body tolerance, low catheter adhesion and radiopacity properties.

Dwg.0/0

TECH UPTX: 20001006

TECHNOLOGY FOCUS - ORGANIC CHEMISTRY - Preferred Monomer: The alkyl **cyanoacrylate** is 2-hexyl **cyanoacrylate**.

Preferred Inhibitors: The monomer component has at least two inhibitors including hydroquinone (50-150 ppm), p-methoxyphenol (50-150 ppm) and **phosphoric acid** (125-375 ppm).

Preferred Fatty Acid Ester: The alkyl esterified fatty acid is selected from alkyl laurate, alkyl palmitate and stearic acid myristate.

Preferred Opacifying Agent: The opacifying agent is selected from platinum, tantalum, titanium, tungsten, barium sulfate and (preferably) gold. The gold is used in the form of a fine powder with a particle size of no more than 7 microns and preferably no more than 1 microns.

TT TT: COMPOSITION ALKYL **CYANOACRYLATE** MONOMER ADHERE SKIN TISSUE.

=> d 1-24 bib hit

L15 ANSWER 1 OF 24 MEDLINE DUPLICATE 1
AN 2002278616 MEDLINE
DN 22000672 PubMed ID: 12006271
TI N-butyl **cyanoacrylate** embolization of cerebral
arteriovenous malformations: results of a prospective,
randomized, multi-center trial.
AU Anonymous
CS The n-BCA Trail Investigators.
SO AJNR. AMERICAN JOURNAL OF NEURORADIOLOGY, (2002 May) 23 (5) 748-55.
Journal code: 8003708. ISSN: 0195-6108.
CY United States
DT (CLINICAL TRIAL)
Journal; Article; (JOURNAL ARTICLE)
(MULTICENTER STUDY)
(RANDOMIZED CONTROLLED TRIAL)
LA English
FS Priority Journals
EM 200206
ED Entered STN: 20020522
Last Updated on STN: 20020612
Entered Medline: 20020611
TI N-butyl **cyanoacrylate** embolization of cerebral
arteriovenous malformations: results of a prospective,
randomized, multi-center trial.
AB BACKGROUND AND PURPOSE: Liquid N-butyl **cyanoacrylate** (n-BCA) use
for the treatment of **arteriovenous malformations** (**AVM**) in the brain has become part of medical practice. However, no
study has led to the Food and Drug Administration's approval of n-BCA for
intravascular use. The purpose of this study was to verify the
effectiveness and safety of an n-BCA/**Tantalum** Powder/Ethiodized
Oil mixture, compared with conventional treatment (Trufill polyvinyl
alcohol [PVA]) for preoperative embolization of cerebral **AVM**.
METHODS: Between October 15, 1996, and March 24, 1999, 104 patients at 13
centers were prospectively randomized to undergo embolization using an
n-BCA/**Tantalum** Powder/Ethiodol mixture or Trufill PVA. The
pre-embolization therapy goals were determined in terms of the number of
pedicles to be embolized and the percent of nidus reduction expected.
Embolization results were evaluated by a central laboratory. Subsequent
surgical resection data were recorded. Safety evaluation data included
recording device complications, procedure complications, and intracranial
events/overall neurologic outcomes, which could be either device-related,
procedure-related, or both. RESULTS: The reduction of **AVM**
dimensions (79.4% in the n-BCA group and 86.9% in the PVA group) and the
mean number of vessels embolized (2.2 in the n-BCA group and 2.1 in the
PVA group) was similar in the two groups. Coils were used more commonly
with PVA embolization (P<.0001). No differences were detected in surgical
resection time, number of patients who required transfusion, volume and
number of transfusion units, or type and volume of fluid replacement.
Glasgow Outcome Scale scores were not significantly different between the
two groups before treatment, after embolization, or after resection. Two
of 42 patients who underwent resection and had been treated with n-BCA
experienced post-resection hematoma, compared with eight of 45 patients
who underwent resection and had been treated with PVA (P<.05). CONCLUSION:
This prospective, randomized trial showed that n-BCA is equivalent to PVA
as a preoperative embolic agent for treatment of cerebral **AVM** as
determined by percent of nidus reduction and number of feeding pedicles
embolized.

CT Check Tags: Comparative Study; Female; Human; Male
 Adult
 *Embolization, Therapeutic
 Enbucrilate: AE, adverse effects
 *Enbucrilate: TU, therapeutic use
 Hematoma: CI, chemically induced
 Hematoma: ET, etiology
 *Intracranial Arteriovenous Malformations: TH, therapy
 Polyvinyl Alcohol: AE, adverse effects
 *Polyvinyl Alcohol: TU, therapeutic use
 Postoperative Complications
 Preoperative Care
 Prospective Studies
 Safety
 Single-Blind Method
 Treatment Outcome

L15 ANSWER 2 OF 24 MEDLINE DUPLICATE 2
 AN 2002114282 IN-PROCESS
 DN 21834597 PubMed ID: 11846039
 TI Combined therapy of cerebral **arteriovenous malformations**
 : histological differences between a non-adhesive liquid embolic agent and
 n-butyl 2-**cyanoacrylate** (NBCA).
 AU Duffner F; Ritz R; Bornemann A; Freudenstein D; Wiendl H; Siekmann R
 CS Department of Neurosurgery, University Hospital, Eberhard-Karls
 University, Tübingen, Germany.. frank.duffner@med.uni-tuebingen.de
 SO CLINICAL NEUROPATHOLOGY, (2002 Jan-Feb) 21 (1) 13-7.
 Journal code: 8214420. ISSN: 0722-5091.
 CY Germany: Germany, Federal Republic of
 DT Journal; Article; (JOURNAL ARTICLE)
 LA English
 FS IN-PROCESS; NONINDEXED; Priority Journals
 ED Entered STN: 20020216
 Last Updated on STN: 20020216
 TI Combined therapy of cerebral **arteriovenous malformations**
 : histological differences between a non-adhesive liquid embolic agent and
 n-butyl 2-**cyanoacrylate** (NBCA).
 AB OBJECTIVE: Based on 2 casuistics, the intraoperative qualities of a new,
 non-adhesive liquid embolic agent (Onyx, Micro Therapeutics. Inc., Irvine,
 CA, USA) are to be compared to those of n-butyl 2-**cyanoacrylate**
 (NBCA) with regard to the histopathological results after preoperative
 embolization of a cerebral **arteriovenous malformation**
 (AVM). PATIENTS AND METHODS: In a case example, the
 intraoperative quality of the nidus after embolization of a
 parieto-occipital **AVM** with Onyx--a new, non-adhesive liquid
 embolic agent--consisting of ethylene-vinyl alcohol copolymer (EVOH),
 dimethyl sulfoxide (DMSO) and **tantalum**, is described. In the
 second patient, embolization of a frontal high-flow **AVM** was
 performed with NBCA. Both patients underwent surgery with complete
 resection of the **AVM**. RESULTS: From a neurosurgical point of
 view, Onyx is suitable for preoperative embolization of **AVMs**,
 because the nidus intraoperatively remains elastic and formable and can be
 dissected from the surrounding brain tissue quite well by microsurgical
 technique. Inflammatory reactions can be found mainly in the lumina of the
 vessels. CONCLUSIONS: Onyx promises to be an embolic agent well suitable
 for subsequent neurosurgical resection. Further studies considering
 various intervals of time between embolization and resection as well as
 histopathological and electron microscopical examinations are necessary
 for evaluation of our first experience with this new embolization agent.

L15 ANSWER 3 OF 24 EMBASE COPYRIGHT 2002 ELSEVIER SCI. B.V.
 AN 2001290389 EMBASE
 TI Successful transarterial glue embolisation by wedged technique for a tentorial dural arteriovenous fistula presenting with a conjunctival injection.
 AU Iizuka Y.; Maehara T.; Hishii M.; Miyajima M.; Arai H.
 CS Y. Iizuka, Department of Radiology, Juntendo University, School of Medicine, 2-1-1 Hongo, Bunkyo-ku, Tokyo 113-8421, Japan
 SO Neuroradiology, (2001) 43/8 (677-679).
 Refs: 15
 ISSN: 0028-3940 CODEN: NRDYAB
 CY Germany
 DT Journal; Article
 FS 008 Neurology and Neurosurgery
 014 Radiology
 037 Drug Literature Index
 LA English
 SL English
 AB Many tentorial dural arteriovenous fistulae (TDAVF) present with intracranial haemorrhage. We report a patient who presented with conjunctival injection. Transarterial embolisation of the TDAVF was undertaken with a wedged injection of a low concentration of N-butyl **cyanoacrylate**, arresting the flow next to the proximal segment of the venous outlet. After three sessions, a complete cure was achieved. We present a useful method which has not been reported previously.
 CT Medical Descriptors:
 *artificial embolism
 *brain arteriovenous malformation: DI, diagnosis
 *brain arteriovenous malformation: DT, drug therapy
 *brain arteriovenous malformation: TH, therapy
 conjunctiva disease
 technique
 cerebellum
 dura mater
 brain hemorrhage
 treatment outcome
 human
 male
 case report
 adult
 article
 priority journal
 Drug Descriptors:
 *embucrilate: CB, drug combination
 *embucrilate: DT, drug therapy
 tissue adhesive: DT, drug therapy
 iodinated poppyseed oil: CB, drug combination
 iodinated poppyseed oil: DT, drug therapy
 tantalum: CB, drug combination
 tantalum: DT, drug therapy
 RN (embucrilate) 25154-80-7, 6606-65-1; (iodinated poppyseed oil) 8001-40-9, 8002-46-8, 8006-56-2, 8006-57-3; (**tantalum**) 7440-25-7

 L15 ANSWER 4 OF 24 CAPLUS COPYRIGHT 2002 ACS DUPLICATE 3
 AN 2000:824137 CAPLUS
 DN 134:9343
 TI Methods for treating **arteriovenous malformations** using radioactive compositions

IN Wallace, George
PA Micro Therapeutics, Inc., USA; Greff, Richard J.
SO PCT Int. Appl., 30 pp.
CODEN: PIXXD2
DT Patent
LA English
FAN.CNT 1

	PATENT NO.	KIND	DATE	APPLICATION NO.	DATE
PI	WO 2000069474	A1	20001123	WO 2000-US13245	20000512
	W:	AE, AG, AL, AM, AT, AU, AZ, BA, BB, BG, BR, BY, CA, CH, CN, CR, CU, CZ, DE, DK, DM, DZ, EE, ES, FI, GB, GD, GE, GH, GM, HR, HU, ID, IL, IN, IS, JP, KE, KG, KP, KR, KZ, LC, LK, LR, LS, LT, LU, LV, MA, MD, MG, MK, MN, MW, MX, NO, NZ, PL, PT, RO, RU, SD, SE, SG, SI, SK, SL, TJ, TM, TR, TT, TZ, UA, UG, US, UZ, VN, YU, ZA, ZW, AM, AZ, BY, KG, KZ, MD, RU, TJ, TM			
	RW:	GH, GM, KE, LS, MW, SD, SL, SZ, TZ, UG, ZW, AT, BE, CH, CY, DE, DK, ES, FI, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE, BF, BJ, CF, CG, CI, CM, GA, GN, GW, ML, MR, NE, SN, TD, TG			
	US 6333020	B1	20011225	US 1999-311803	19990513
	EP 1191947	A1	20020403	EP 2000-932420	20000512
	R:	AT, BE, CH, DE, DK, ES, FR, GB, GR, IT, LI, LU, NL, SE, MC, PT, IE, SI, LT, LV, FI, RO			
PRAI	US 1999-311803	A1	19990513		
	WO 2000-US13245	W	20000512		

RE.CNT 3 THERE ARE 3 CITED REFERENCES AVAILABLE FOR THIS RECORD
ALL CITATIONS AVAILABLE IN THE RE FORMAT

TI Methods for treating **arteriovenous malformations** using radioactive compositions

AB Disclosed are methods for treating **arteriovenous malformations** in a mammal by use of a radioactive compn. The methods involve in vivo delivery of radioactive substances in a fluid to one or more vascular sites in the **arteriovenous malformation**. Subsequent solidification of the compn. results in vascular embolization to partially ablate the **arteriovenous malformation** and delivery of a controlled amt. of radiation to further ablate the **arteriovenous malformation** and inhibits its regrowth. In one example, a compn., using a "cold" isotope, is prepd. based on ethylene-vinyl alc. copolymer, micronized **tantalum**, iridium powder, and DMSO. Such compns., when added to saline, provide a solid, coherent ppt. The delivery of such compns. was demonstrated in pigs. A non-radioactive contrast agent may be added to the compns. to enable visualization of the delivery.

ST polymeric radioactive compn **arteriovenous malformation** ablation

IT Blood vessel, disease
(**arteriovenous malformation**; radioactive compns. for treating **arteriovenous malformations**)

IT Solvents
(biocompatible; radioactive compns. for treating **arteriovenous malformations**)

IT Polymers, biological studies
RL: BPR (Biological process); BSU (Biological study, unclassified); FMU (Formation, unclassified); PEP (Physical, engineering or chemical process); THU (Therapeutic use); BIOL (Biological study); FORM (Formation, nonpreparative); PROC (Process); USES (Uses)
(biocompatible; radioactive compns. for treating **arteriovenous malformations**)

IT Imaging agents

(contrast; radioactive compns. for treating **arteriovenous malformations**)

IT Polyesters, biological studies
 RL: THU (Therapeutic use); BIOL (Biological study); USES (Uses)
 (glycolide-based; radioactive compns. for treating **arteriovenous malformations**)

IT Polyesters, biological studies
 RL: THU (Therapeutic use); BIOL (Biological study); USES (Uses)
 (lactide; radioactive compns. for treating **arteriovenous malformations**)

IT Polyethers, biological studies
 RL: THU (Therapeutic use); BIOL (Biological study); USES (Uses)
 (ortho ester group-contg.; radioactive compns. for treating **arteriovenous malformations**)

IT Polyesters, biological studies
 RL: THU (Therapeutic use); BIOL (Biological study); USES (Uses)
 (polyamide-; radioactive compns. for treating **arteriovenous malformations**)

IT Polyethers, biological studies
 RL: THU (Therapeutic use); BIOL (Biological study); USES (Uses)
 (polycarbonate-; radioactive compns. for treating **arteriovenous malformations**)

IT Polyamides, biological studies
 RL: THU (Therapeutic use); BIOL (Biological study); USES (Uses)
 (polyester-; radioactive compns. for treating **arteriovenous malformations**)

IT Polycarbonates, biological studies
 RL: THU (Therapeutic use); BIOL (Biological study); USES (Uses)
 (polyether-; radioactive compns. for treating **arteriovenous malformations**)

IT Drug delivery systems
 Drug targeting
 Hydrogels
 Radiopharmaceuticals
 Radiotherapy
 (radioactive compns. for treating **arteriovenous malformations**)

IT Radionuclides, biological studies
 RL: BPR (Biological process); BSU (Biological study, unclassified); THU (Therapeutic use); BIOL (Biological study); PROC (Process); USES (Uses)
 (radioactive compns. for treating **arteriovenous malformations**)

IT Collagens, biological studies
 Fibrins
 Gelatins, biological studies
 Polyamides, biological studies
 Polyanhydrides
 Polycarbonates, biological studies
 Polyesters, biological studies
 Polyketones
 Polyoxymethylenes, biological studies
 Polyoxymethylenes, biological studies
 Polyphosphazenes
 Polyurethanes, biological studies
 RL: THU (Therapeutic use); BIOL (Biological study); USES (Uses)
 (radioactive compns. for treating **arteriovenous malformations**)

IT 440-58-4, Iodamide 1225-20-3, Iothalamate sodium 1314-61-0,
 Tantalum oxide 6284-40-8, Meglumine 7440-06-4,

Platinum, biological studies 7440-25-7, **Tantalum**, biological studies 7440-33-7, **Tungsten**, biological studies 7440-57-5, **Gold**, biological studies 7727-43-7, **Barium sulfate** 31112-62-6, **Metrizamide** 60166-93-0, **Iopamidol**

RL: THU (Therapeutic use); BIOL (Biological study); USES (Uses)
(contrast agent; radioactive compns. for treating **arteriovenous malformations**)

IT 25154-80-7, Poly(butyl **cyanoacrylate**)
RL: BPR (Biological process); BSU (Biological study, unclassified); FMU (Formation, unclassified); PEP (Physical, engineering or chemical process); THU (Therapeutic use); BIOL (Biological study); FORM (Formation, nonpreparative); PROC (Process); USES (Uses)

(radioactive compns. for treating **arteriovenous malformations**)

IT 6606-65-1
RL: BPR (Biological process); BSU (Biological study, unclassified); PEP (Physical, engineering or chemical process); RCT (Reactant); THU (Therapeutic use); BIOL (Biological study); PROC (Process); RACT (Reactant or reagent); USES (Uses)

(radioactive compns. for treating **arteriovenous malformations**)

IT 67-68-5, **Dmsol**, biological studies 7439-88-5, **Iridium**, biological studies 25067-34-9, **Ethylene-vinyl alcohol copolymer**

RL: BPR (Biological process); BSU (Biological study, unclassified); PEP (Physical, engineering or chemical process); THU (Therapeutic use); BIOL (Biological study); PROC (Process); USES (Uses)

(radioactive compns. for treating **arteriovenous malformations**)

IT 1398-61-4, **Chitin** 9003-20-7, **Polyvinyl acetate** 9003-39-8, **Polyvinylpyrrolidone** 9004-34-6, **Hydroxycellulose**, biological studies 9004-35-7, **Cellulose acetate** 9004-36-8, **Cellulose acetate butyrate** 9004-70-0, **Nitrocellulose** 9012-76-4, **Chitosan** 10043-49-9, **Gold** 198, biological studies 10045-97-3, **Cesium** 137, biological studies 10098-91-6, **Yttrium** 90, biological studies 10098-97-2, **Strontium** 90, biological studies 10198-40-0, **Cobalt** 60, biological studies 13981-50-5, **Cobalt** 57, biological studies 13981-99-2, **Nickel** 57, biological studies 13982-25-7, **Cobalt** 55, biological studies 14093-03-9, **Cobalt** 56, biological studies 14093-04-0, **Iron** 52, biological studies 14119-09-6, **Gallium** 67, biological studies 14158-31-7, **Iodine** 125, biological studies 14391-22-1, **Thulium** 167, biological studies 14596-37-3, **Phosphorus** 32, biological studies 14681-59-5, **Iron** 55, biological studies 14687-25-3, **Lead** 203, biological studies 14694-69-0, **Iridium** 192, biological studies 14809-46-2, **Selenium** 72, biological studies 14833-23-9, **Zinc** 62, biological studies 14967-68-1, **Palladium** 103, biological studies 15047-05-9, **Cesium** 129, biological studies 15064-65-0, **Thallium** 201, biological studies 15128-03-7, **Copper** 61, biological studies 15422-57-8, **Selenium** 73, biological studies 15715-08-9, **Iodine** 123, biological studies 15741-33-0, **Manganese** 57, biological studies 15750-15-9, **Indium** 111, biological studies 15755-33-6, **Arsenic** 72, biological studies 15765-39-6, **Bromine** 77, biological studies 15776-19-9, **Bismuth** 206, biological studies 18268-34-3, **Rubidium** 81, biological studies 24980-41-4, **Polycaprolactone** 25014-41-9, **Polyacrylonitrile** 25248-42-4, **Polycaprolactone** 25300-64-5, **Styrene-maleic acid copolymer** 25322-68-3, **Polyethylene glycol** 26009-03-0, **Polyglycolide** 26023-30-3, **Poly[oxy(1-methyl-2-oxo-1,2-ethanediyl)]** 26063-00-3, **Polyhydroxybutyrate** 26202-08-4, **Polyglycolide** 26680-10-4, **Poly lactide** 26744-04-7 31621-87-1, **Polydioxanone** 78644-42-5, **Poly(malic acid)** 102190-94-3, **Polyhydroxyvaleric acid**

RL: THU (Therapeutic use); BIOL (Biological study); USES (Uses)
(radioactive compns. for treating **arteriovenous malformations**)

IT 64-17-5, Ethanol, biological studies 67-64-1, Acetone, biological studies 97-64-3, Ethyl lactate

RL: THU (Therapeutic use); BIOL (Biological study); USES (Uses)
(solvent; radioactive compns. for treating **arteriovenous malformations**)

L15 ANSWER 5 OF 24 WPIDS (C) 2002 THOMSON DERWENT

AN 2001-070908 [08] WPIDS

DNN N2001-053688 DNC C2001-019754

TI Intracorporeal space filling device for treating e.g. patient's blood vessels, includes material which is transmutable from non-rigid state to rigid state.

DC A14 A17 A28 A96 D22 P32

IN MARKS, M P; ROSS, M

PA (SETH-N) SETHEL INTERVENTIONAL INC

CYC 93

PI WO 2000072781 A2 20001207 (200108)* EN 40p

RW: AT BE CH CY DE DK EA ES FI FR GB GH GM GR IE IT KE LS LU MC MW MZ
NL OA PT SD SE SL SZ TZ UG ZW

W: AE AG AL AM AT AU AZ BA BB BG BR BY CA CH CN CR CU CZ DE DK DM DZ
EE ES FI GB GD GE GH GM HR HU ID IL IN IS JP KE KG KP KR KZ LC LK
LR LS LT LU LV MA MD MG MK MN MW MX MZ NO NZ PL PT RO RU SD SE SG
SI SK SL TJ TM TR TT TZ UA UG UZ VN YU ZA ZW

AU 2000051799 A 20001218 (200118)

EP 1200012 A2 20020502 (200236) EN

R: AL AT BE CH CY DE DK ES FI FR GB GR IE IT LI LT LU LV MC MK NL PT
RO SE SI

ADT WO 2000072781 A2 WO 2000-US15445 20000602; AU 2000051799 A AU 2000-51799
20000602; EP 1200012 A2 EP 2000-936492 20000602, WO 2000-US15445 20000602

FDT AU 2000051799 A Based on WO 200072781; EP 1200012 A2 Based on WO 200072781

PRAI US 1999-324987 19990602

AB WO 200072781 A UPAB: 20010207

NOVELTY - An intracorporeal space filling device comprises a transmutable material disposed within an inner lumen (15) of an elongate tubular shell. The inner lumen is fluidly connected with first and second ports which are disposed at first and second ends of the shell, respectively. The material is transmutable from a non-rigid to a rigid state within the patient's body.

DETAILED DESCRIPTION - INDEPENDENT CLAIMS are also included for:

(a) a method of **occluding** an intracorporeal void comprising positioning the distal end of a microcatheter so that a distal port in the distal end will be directed to the cavity of the intracorporeal void, advancing an intracorporeal space filling device (10), and transmuting the transmutable material (16); and

(b) a detachment mechanism for detaching an intracorporeal space filling device from a delivery system comprising a degradable polymer link securing the intracorporeal space filling device to the delivery system, and a heating element disposed in thermal contact with the degradable polymer link.

USE - The device is used to treat a patient's blood vessels, intracorporeal conduits or other portions of the patient's body. It can also be used to treat intracranial aneurysms, arteriovenous fistulas, and other abnormalities within the cerebral vasculature.

ADVANTAGE - The transmutable material provides a space filling device to be soft and flexible at the time of deployment into an intracorporeal cavity and is incompressible after converted to a rigid state. The device

conforms to the morphology of intracorporeal cavities and transmutes to rigid mass upon activation or hardening of the transmutable material. The device is resistant to compression and reforming due to vascular or other types of pressures within the patient's body.

DESCRIPTION OF DRAWING(S) - The figure is a longitudinal sectional view of an intracorporeal space filling device.

Intracorporeal space filling device 10

Shell 11

Lumen 15

Transmutable material 16

Dwg.1/26

TECH

UPTX: 20010207

TECHNOLOGY FOCUS - POLYMERS - Preferred Materials: The outer wall of the shell (11) is a polymer material of polyurethane, polyethylene, nylon, polyimide, polyamide, polytetrafluoroethylene, polyester or polypropylene. The transmutable material is methacrylate compounds, linear polyester, silicone, **cyanoacrylates**, polyisocyanate, ultra violet curable acrylates, moisture cure silicones, dimethyl sulfoxide, thioisocyanate aldehyde, isocyanate, divinyl compounds, epoxide acrylates, succinimidyl azido salicylate, succinimidyl azidobenzoate, succinimidyl dithio acetate, azidoiodobenzene, fluoronitrophenylazide, salicylate azides or benzophenonemaleimide.

TECHNOLOGY FOCUS - METALLURGY - Preferred Material: The shell may be a metal of stainless steel, nickel titanium alloy, **gold**, **platinum**, **tantalum** or palladium.

TECHNOLOGY FOCUS - INSTRUMENTATION AND TESTING - Preferred Components: The shell has an outside surface, which is self adhering in a fluid field to create attachment from contact points upon activation. The shell has apertures to expose portions of the transmutable material upon transmutation. An elongated longitudinal mechanism is secured to the device coextensive with the shell. A helical coil is disposed to, and coextensive with the shell. A delivery system, the first end of the shell, and bead(s) containing transmutable material are detachably secured to the detachment mechanism. Each bead is connected to adjacent bead by a flexible mechanism and configured to produce a linear array of the beads. The flexible mechanism is a portion of an elongated mechanism disposed along axis of the device.

Preferred Method: A blocking balloon, which is adjacent to the intracorporeal void and the distal end of the microcatheter, is deployed before advancing the device into the void. The polymer link is degraded by a chain cleavage reaction. The heating element is configured to be heated by an electric current.

L15 ANSWER 6 OF 24 WPIDS (C) 2002 THOMSON DERWENT

AN 2000-602071 [57] WPIDS

CR 2000-223969 [18]; 2000-572323 [53]; 2002-214894 [26]

DNN N2000-445484 DNC C2000-180201

TI Magnetic embolic agent for treatment of vascular defects, particularly an aneurysm or atriovenous malformation, comprises a polymer, solvent, adhesive and magnetic particles.

DC A18 A25 A96 B04 D22 G03 P34 S05

IN GARIBALDI, J M; HASTINGS, R N; HOGG, B J; REN, B

PA (STER-N) STEREOTAXIS INC

CYC 91

PI WO 2000054832 A1 20000921 (200057)* EN 56p

RW: AT BE CH CY DE DK EA ES FI FR GB GH GM GR IE IT KE LS LU MC MW NL

OA PT SD SE SL SZ TZ UG ZW

W: AE AL AM AT AU AZ BA BB BG BR BY CA CH CN CR CU CZ DE DK DM EE ES

FI GB GD GE GH GM HR HU ID IL IN IS JP KE KG KP KR KZ LC LK LR LS
LT LU LV MA MD MG MK MN MW MX NO NZ PL PT RO RU SD SE SG SI SK SL
TJ TM TR TT TZ UA UG UZ VN YU ZA ZW

AU 2000038983 A 20001004 (200101)

US 6296604 B1 20011002 (200160)

EP 1169081 A1 20020109 (200205) EN

R: AL AT BE CH CY DE DK ES FI FR GB GR IE IT LI LT LU LV MC MK NL PT
RO SE SI

ADT WO 2000054832 A1 WO 2000-US7222 20000316; AU 2000038983 A AU 2000-38983
20000316; US 6296604 B1 CIP of US 1999-271118 19990317, US 1999-430200
19991029; EP 1169081 A1 EP 2000-918116 20000316, WO 2000-US7222 20000316

FDT AU 2000038983 A Based on WO 200054832; EP 1169081 A1 Based on WO 200054832

PRAI US 1999-430200 19991029; US 1999-271118 19990317

AB WO 200054832 A UPAB: 20020429

NOVELTY - A magnetic embolic agent (MEA) for magnetic placement in a
vascular defect (VD) to form an embolus in it to **occlude** it, is
new.

DETAILED DESCRIPTION - A novel MEA for magnetic placement in a VD to
form an embolus in it to **occlude** the defect, comprises:

- (a) 5 - 50 wt. % biocompatible polymer (BP);
- (b) 40 - 90 wt. % biocompatible solvent (BS) capable of solubilizing
the BP;
- (c) 1 - 20 wt. % adhesive (AA); and
- (d) 5 - 50 wt. % magnetic particles (MP) responsive to a magnetic
field.

INDEPENDENT CLAIMS are also included for the following:

(1) treating a VD comprising magnetically positioning a MEA at the
site of the VD to form an embolus at the site of the VD; the MEA being
chemically reactive to lose some of its magnetic responsiveness after
formation of the embolus;

(2) treating a VD comprising:
(a) providing a real time digital image of the operating site;
(b) applying a magnetic field at the operating site with an external
magnet;

(c) ejecting a magnetic embolic material from a catheter under the
influence of the applied magnetic field to form an embolus that
occludes the defect; and

(d) stopping the ejection of embolic material when the image
indicates that the defect is **occluded**;

(3) a biocompatible magnetic mixture (BMM) which can be ejected
through a standard microcatheter (MC) and held magnetically within a VD,
comprising:

- (a) 5-50 wt.% MP;
- (b) 5-50wt.% BP;
- (c) 40-90wt.% solvent which dissolves the BP; the MP and BP being
homogeneously distributed throughout the BS which precipitates from the
solution as the BMM is deposited into the VD and held by an externally
applied magnetic field;

(4) a BMM which can be delivered through a MC and held magnetically
within a VD composed of:

- (a) 5-50wt.% MP;
- (b) 5-50wt.% BP;
- (c) 1-20wt.% AA;
- (d) 40-90wt.% BS which dissolves the BP; the AA improving the
cohesiveness and tissue attachment characteristics of the BMM while the
MP, AA and BP are homogeneously distributed throughout the BS and
precipitates from the solution as the BMM is deposited into the VD and
held by an externally applied magnetic field;

(5) a BMM which can be delivered through a MC and held magnetically

within a VD composed of:

- (a) 5-50wt.% MP; and
- (b) 5-50wt.% AA; the MP isomogeneously distributed within the BMM and the BMM deposited into the VD and held by an externally applied magnetic field;

(6) a BMM can be delivered through a MC and held magnetically within a VD composed of:

- (a) 5-50wt.% MP;
- (b) 1-50wt.% dispersion agent;
- (c) 10-95wt.% AA; where the dispersion agent improves the homogeneity of the MP within the BMM and the BMM can be deposited into the VD and held by an externally applied magnetic field;

(7) a BMM which can be delivered through a MC and held magnetically within a VD composed of:

- (a) 5-90wt.% coated MP;
- (b) 10-90wt.% BS; whereby the MP is coated with an agent which is dissolved in the BS to form a homogeneous BMM which can be deposited into the VD and held by an externally applied magnetic field;

(8) a MEA for magnetic placement in a VD with increased X-ray opacification to form an embolus in the defect to occlude the defect, the agent comprising:

- (a) 4-70wt.% BP;
- (b) 10-80wt.% BS capable of solubilizing the BP;
- (c) 10-50wt.% MP responsive to a magnetic field; and
- (d) 10-50wt.% X-ray opaque MP responsive to a magnetic field;

(9) a MEA for magnetic placement in a VD using diluted solvents to form an embolus in the defect to occlude the defect, the agent comprising:

- (a) 4-70wt.% biocompatible reactive polymer;
- (b) 10-80wt.% BS diluted in water capable of solubilizing the BP;
- (c) 0-50wt.% BP; and
- (d) 10-50wt.% MP responsive to a magnetic field;

(10) a MEA for magnetic placement in a VD using diluted solvents with increased X-ray opacification to form an embolus in the defect to occlude the defect, the agent comprising:

- (a) 4-70wt.% BP;
- (b) 10-80wt.% BS diluted in water capable of solubilizing the BP;
- (c) 0-50wt.% adhesive;
- (d) 10-50wt.% MP responsive to a magnetic field; and
- (e) 10-50wt.% X-ray opaque MP responsive to a magnetic field;

(11) an embolic agent for delivery into a VD to form an embolus in the defect to occlude the defect, the agent comprising:

- (a) 4-80wt.% BP;
- (b) 30-95wt.% BS capable of solubilizing the BP;
- (c) 1-70wt.% adhesive; and
- (d) an X-ray opaque material is added to enhance the visibility under fluoroscopy;

(12) an embolic agent for delivery into a VD to form an embolus in the defect to occlude the defect, the agent comprising:

- (a) 4-80wt.% biocompatible reactive polymer;
- (b) 10-90wt.% BS diluted in water capable of solubilizing the BP;
- (c) 0-80wt.% BP; and
- (d) an X-ray opaque material to enhance visibility under fluoroscopy;

(13) a MEA for magnetic placement in a VD to form an embolus in it to occlude the defect, the agent comprising:

- (a) 10-90wt.% biocompatible reactive polymer;
- (b) 10-80wt.% MP responsive to a magnetic field; and
- (c) 10-80wt.% X-ray opaque MP responsive to a magnetic field;

(14) a 2-part MEA for magnetic placement in a VD to form an embolus in it to occlude the defect, the agent comprising:

- (a) a first part comprising:
 - (i) 10-90wt.% biocompatible reactive polymer;
 - (ii) 10-80wt.% MP responsive to a magnetic field; and
 - (iii) 10-80wt.% X-ray opaque MP responsive to a magnetic field; and
- (b) a second part comprising 10-90wt.% BP catalyst;
- (15) treating a VD comprising:
 - (a) introducing a flowable first magnetic composition into the VD under the guidance of an externally applied magnetic field;
 - (b) introducing a flowable second magnetic composition into the VD under the guidance of an externally applied magnetic field, the second magnetic composition when mixed with the first magnetic composition forming a non-flowable material; and
 - (c) mixing the first and second magnetic compositions in the VD by varying the externally applied magnetic field to form an occlusion in the VD;
- (16) retarding the hardening of an embolic material injected into a VD comprising injecting a biocompatible liquid with a high surface tension prior to injecting the embolic material to create a clean barrier between bodily fluid and embolic material;
- (17) a magnetic liquid embolic agent responsive to an externally applied magnetic field to flow into a VD and harden to occlude the VD, the embolic agent comprising a BP, a BS, 25-40wt.% magnetite, and 15-25wt.% gold plated nickel.

ACTIVITY - Vasotropic. No biological data is given.

MECHANISM OF ACTION - None given.

USE - The new compositions and methods can be used for treating vascular defects such as aneurysms and atriovenous malformations.

ADVANTAGE - The new compositions can be delivered intravascularly through a catheter and can be guided into and held in place in the vascular defect with an applied magnetic field. Defects of all shapes and at all locations can be treated equally by simply adjusting the magnetic force direction.

Dwg.21/23

TECH

UPTX: 20001109

TECHNOLOGY FOCUS - POLYMERS - Preferred Agent: The BP includes polyurethane, hydrogel, ethylene vinyl alcohol (EVOH), polynethylmethacrylate, cellulose acetate, polyvinyl alcohol, prolamines, ethyl cellulose, polyvinyl acetate, polyvinyl butyrate, polyvinyl alcohol, hydrogels, polyvinyl pyrrolidone, or mussel adhesive protein. The coating for the BMM includes at least one of cellulose acetate, polymerized methylmethacrylate, polyethylene vinyl alcohol, or polyvinyl acetate. In the magnetic embolic agents, the biocompatible polymer may comprise prolamine, and the diluting solvent may comprise ethanol.

TECHNOLOGY FOCUS - ORGANIC CHEMISTRY - Preferred Agent: The BS is preferably dimethylsulphoxide (DMSO), acetone, ethyl acetate or ethanol. The AA is preferably **cyanoacrylate**, methylmethacrylate or fibrin. The ratio of the BP to the biocompatible adhesive may be 40:1 to 1:1. The dispersion agent for the BMM may be at least one of a stearate, seed oils or esters.

TECHNOLOGY FOCUS - INORGANIC CHEMISTRY - Preferred Agent: The MP may be 1-10 micrometers and are preferably iron, coated iron, carbonyl-iron, or iron carbon composite particles.

L15 ANSWER 7 OF 24 MEDLINE

DUPLICATE 4

AN 2000248775 MEDLINE

DN 20248775 PubMed ID: 10789913

TI Pulmonary emboli following therapeutic embolization of cerebral

arteriovenous malformations in children.

AU Kjellin I B; Boechat M I; Vinuela F; Westra S J; Duckwiler G R
CS Department of Radiological Sciences, UCLA School of Medicine, Center of
the Health Sciences, Los Angeles, CA 90095-1721, USA.
SO PEDIATRIC RADIOLOGY, (2000 Apr) 30 (4) 279-83.
Journal code: 0365332. ISSN: 0301-0449.
CY GERMANY: Germany, Federal Republic of
DT Journal; Article; (JOURNAL ARTICLE)
LA English
FS Priority Journals
EM 200006
ED Entered STN: 20000622
Last Updated on STN: 20000622
Entered Medline: 20000609
TI Pulmonary emboli following therapeutic embolization of cerebral
arteriovenous malformations in children.
AB BACKGROUND: Reports of the complicating side effect of pulmonary embolism
(PE) following endovascular therapy of cerebral **arteriovenous**
malformations (AVM) in children have been limited in
number. Details of its occurrence are yet to be fully elucidated.
OBJECTIVE: The hypothesis is that inadvertent pulmonary migration of
embolic material is common and may go unrecognized. MATERIALS AND METHODS:
Forty-seven patients (ages 1 day to 16 years and 11 months) underwent
embolization of a cerebral **AVM** with at least one material (
cyanoacrylate, **platinum** coils, detachable balloons,
polyvinyl alcohol particles). The medical records and chest radiographs
were reviewed retrospectively. Chest radiographs were available in 34
patients. The radiographs were analyzed for the presence or absence of
foreign material in the lungs. RESULTS: The chest radiographs in 12
patients (35%) showed pulmonary deposits of embolic material;
cyanoacrylate in 10 patients and **platinum** coils in 2.
Two of the patients with **cyanoacrylate** deposits in the lungs
developed respiratory distress that required endotracheal intubation. The
patients gradually improved after a time period of 7-10 days with
conservative treatment. CONCLUSION: PE is not an uncommon complication in
children undergoing embolization of brain **AVM**. Although usually
asymptomatic, PE may cause severe symptoms.
CT Check Tags: Comparative Study; Female; Human; Male
Adolescence
Child
Child, Preschool
Electrocardiography
*Embolization, Therapeutic: AE, adverse effects
Embolization, Therapeutic: MT, methods
Enbucrilate
Infant
Infant, Newborn
*Intracranial Arteriovenous Malformations: TH, therapy
Platinum
Polyvinyl Alcohol
Pulmonary Embolism: DI, diagnosis
*Pulmonary Embolism: ET, etiology
Pulmonary Embolism: RA, radiography
Radiography, Thoracic
Tomography, X-Ray Computed
RN 6606-65-1 (Enbucrilate); 7440-06-4 (Platinum); 9002-89-5
(Polyvinyl Alcohol)

AN 1999063588 MEDLINE
 DN 99063588 PubMed ID: 9848842
 TI Embolization of neurosurgical lesions involving the ophthalmic artery.
 AU Lefkowitz M; Giannotta S L; Hieshima G; Higashida R; Halbach V; Dowd C; Teitelbaum G P
 CS Department of Neurological Surgery, University of Southern California School of Medicine, Los Angeles, USA.
 SO NEUROSURGERY, (1998 Dec) 43 (6) 1298-303.
 Journal code: 7802914. ISSN: 0148-396X.
 CY United States
 DT Journal; Article; (JOURNAL ARTICLE)
 LA English
 FS Priority Journals
 EM 199903
 ED Entered STN: 19990326
 Last Updated on STN: 19990326
 Entered Medline: 19990312
 AB OBJECTIVE: A number of anteriorly located cranial base and extracranial lesions receive their vascular supply wholly or in part from the ophthalmic artery, and embolization of the ophthalmic artery can be helpful in the management of these lesions, either as the primary treatment or as an adjunct to surgery. We present situations in which the embolization of lesions involving the ophthalmic artery was performed to effect a partial or total cure of the lesion. METHODS: Twelve patients underwent a total of 15 embolization attempts on lesions involving the ophthalmic artery. Four patients had **arteriovenous malformations** of the orbit, four had dural arteriovenous fistulae, two had orbital meningiomas, one had a planum sphenoidale meningioma, and one had a juvenile nasal angiofibroma. In each case, a Tracker No. 18 microcatheter (Target Therapeutics, Inc., Fremont, CA) was navigated into the ophthalmic artery using a steerable guidewire and digital road mapping. Embolic agents included polyvinyl alcohol particles ranging from 350 to 1500 microm in diameter, 2-mm **platinum** microcoils, and n-butyl-**cyanoacrylate**. In 12 of 15 cases, lidocaine and amytal provocation tests were conducted before any attempt at embolization to assess the role of the ophthalmic artery in vision. RESULTS: Embolization was successfully performed in the 14 situations in which it was attempted. Positive results of two lidocaine/amytal tests were noted. In one case, embolization was not attempted. In the other case, a larger caliber embolic agent (2-mm **platinum** coils) was used. A single transient decrease in visual acuity lasting 4 days was the only embolization-related complication. CONCLUSION: Proper case selection, judicious use of embolic agents, and use of provocative testing can result in safe embolization of lesions supplied by the ophthalmic artery.
 CT Check Tags: Case Report; Female; Human; Male
 Adolescence
 Adult
 Aged
 Amobarbital: DU, diagnostic use
 Angiofibroma: BS, blood supply
 Angiofibroma: SU, surgery
 *Angiofibroma: TH, therapy
 Arteriovenous Fistula: SU, surgery
 Arteriovenous Fistula: TH, therapy
Arteriovenous Malformations: SU, surgery
 ***Arteriovenous Malformations: TH, therapy**
 Combined Modality Therapy
 *Dura Mater: BS, blood supply
 *Embolization, Therapeutic

Embolization, Therapeutic: AE, adverse effects
 Embolization, Therapeutic: IS, instrumentation
 Embolization, Therapeutic: MT, methods
 Enbucrilate: TU, therapeutic use
 Infant
 Lidocaine: DU, diagnostic use
 Meningeal Neoplasms: BS, blood supply
 Meningeal Neoplasms: SU, surgery
 *Meningeal Neoplasms: TH, therapy
 Meningioma: BS, blood supply
 Meningioma: SU, surgery
 *Meningioma: TH, therapy
 Middle Age
 Nose Neoplasms: BS, blood supply
 Nose Neoplasms: SU, surgery
 *Nose Neoplasms: TH, therapy
 *Ophthalmic Artery
 Ophthalmic Artery: AH, anatomy & histology
 Orbital Neoplasms: BS, blood supply
 Orbital Neoplasms: SU, surgery
 *Orbital Neoplasms: TH, therapy
 Particle Size
 Polyvinyl Alcohol: TU, therapeutic use
 Prostheses and Implants
 Retinal Artery: AH, anatomy & histology
 Treatment Outcome
 Vision Disorders: ET, etiology
 Vision Disorders: PC, prevention & control

L15 ANSWER 9 OF 24 EMBASE COPYRIGHT 2002 ELSEVIER SCI. B.V.
 AN 1998374545 EMBASE
 TI [Anaesthesiological management of patients with arteriovenous malformations (**AVMs**) undergoing neuroradiological intervention].
 ANASTHESIOLOGISCHES MANAGEMENT BEI PATIENTEN MIT ARTERIOVENOSEN MALFORMATIONEN (**AVM**) IN DER INTERVENTIONELLEN NEURORADIOLOGIE.
 AU Jaeger K.; Ruschulte H.; Heine J.; Leuwer M.; Piepenbrock S.
 CS Dr. K. Jaeger, Zentrum Anesthesiologie, Medizinische Hochschule Hannover, Carl-Neuberg-Strasse 1, D-30625 Hannover, Germany
 SO Anesthesiologie und Intensivmedizin, (1998) 39/10 (501-504).
 Refs: 21
 ISSN: 0170-5334 CODEN: ANIMD2
 CY Germany
 DT Journal; Article
 FS 008 Neurology and Neurosurgery
 024 Anesthesiology
 LA German
 SL English; German
 TI [Anaesthesiological management of patients with arteriovenous malformations (**AVMs**) undergoing neuroradiological intervention].
 ANASTHESIOLOGISCHES MANAGEMENT BEI PATIENTEN MIT ARTERIOVENOSEN MALFORMATIONEN (**AVM**) IN DER INTERVENTIONELLEN NEURORADIOLOGIE.
 AB Intracranial **arteriovenous malformations (AVM)**
) present with neurological symptoms like headache, focal and general convulsions or disordered vigilance mainly caused by bleeding or infarction. The treatment of **AVMs** consists of neurosurgery, radiosurgery or interventional neuroradiology: Vessels feeding or draining an **AVM** convolute are an **occluded** selectively by N-butyl **cyanoacrylate** or **platinum** coils, probably neccessitating several interventional sessions. Neuroradiological

treatment of **arteriovenous malformations** has been remarkably improved over recent years. Endovascular embolisation can be performed under sedation or general anaesthesia. With respect to the delicate anatomic and pathophysiological condition of **AVMs**, appropriate periinterventional anaesthesiological monitoring and treatment have to be chosen. Intracranial haemodynamics and brain metabolism may not be irritated by drugs and anaesthesia management: For early neurological assessment patients should be wide awake once neuroradiological procedures are finished. Central nervous functions should be monitored postoperatively in an intermediate care or intensive care unit. Basically, principles of neurosurgical anaesthesia can be transferred to anaesthesia management of patients undergoing neuroradiological procedures.

CT Medical Descriptors:

***arteriovenous malformation**
 *neuroradiology
 *anesthesiology
 neurologic disease: DI, diagnosis
 preoperative care
 artificial embolism
 postoperative complication
 neurologic examination
 intensive care
 patient care
 sedation
 general anesthesia
 human
 article

L15 ANSWER 10 OF 24 EMBASE COPYRIGHT 2002 ELSEVIER SCI. B.V.

AN 1999029567 EMBASE

TI Influence of temperature on embolisation with **cyanoacrylate**.

AU Bracard S.; Macho-Fernandez J.M.; Wang X.; Anxionnat R.; Picard L.

CS Prof. L. Picard, University Hospital, 1, Rue Foiler C.O. n. 34, 54035 Nancy Cedex, France. s.bracard@chu-nancy.fr

SO Interventional Neuroradiology, (1998) 4/4 (301-305).

Refs: 15

ISSN: 1123-9344 CODEN: INEUF5

CY Italy

DT Journal; Article

FS 008 Neurology and Neurosurgery

014 Radiology

LA English

SL English

TI Influence of temperature on embolisation with **cyanoacrylate**.

AB We evaluated the influence of temperature on the viscosity of mixtures with different histoacryl/lipiodol concentrations and on injection control, to test the radiological visualization of these mixtures. A viscosimeter was used to measure the viscosity of different histoacryl and lipiodol combinations at various temperatures. After introduction of these blends into the polyethylene tubes, their radiological densities were evaluated by means of CT and DSA. Viscosity was found to be directly proportional to the percentage of lipiodol and inversely proportional to the temperature. By digital subtraction, the mixtures were still visible when the percentage of histoacryl reached 90%. Warming histoacryl and lipiodol mixtures to a temperature that is close to 40.degree.C decreases the mixture's viscosity significantly and makes the injection easier to manage. **Tantalum** and tungsten powders do not necessarily have to be added to visualize mixtures containing less than 90% histoacryl.

CT Medical Descriptors:

*artificial embolism
*brain artery aneurysm: SU, surgery
*brain arteriovenous malformation: SU, surgery

viscosity
computer assisted tomography
viscometry
digital subtraction angiography
temperature sensitivity
article

Drug Descriptors:

*cyanoacrylate

RN (cyanoacrylate) 15802-18-3

L15 ANSWER 11 OF 24 MEDLINE DUPLICATE 6

AN 97022702 MEDLINE

DN 97022702 PubMed ID: 8869062

TI Endovascular treatment of experimental aneurysms with liquid polymers: the protective potential of stents.

AU Szikora I; Guterman L R; Standard S C; Wakhloo A K; Hopkins L N

CS Department of Neurosurgery, State University of New York at Buffalo, USA.

SO NEUROSURGERY, (1996 Feb) 38 (2) 339-47.

Journal code: 7802914. ISSN: 0148-396X.

CY United States

DT Journal; Article; (JOURNAL ARTICLE)

LA English

FS Priority Journals

EM 199701

ED Entered STN: 19970128

Last Updated on STN: 19970128

Entered Medline: 19970114

AB Liquid polymers have previously been used to treat experimental and human aneurysms. However, the delivery of a liquid embolic material into the cerebral circulation involves a high risk of irreversible vessel **occlusion** and stroke. To evaluate methods for the safe and effective treatment of experimental aneurysms with liquid polymer injection, we tested four different techniques to deliver cellulose acetate polymer (CAP) or N-hexyl-**cyanoacrylate** into canine side-wall carotid artery aneurysms. The animals were observed for 1 to 10 weeks after treatment. Two aneurysms were treated without protection of the distal circulation, one with CAP and another with N-hexyl-**cyanoacrylate**. In four cases, an angioplasty balloon was inflated within the parent artery during endosaccular injection of CAP. In two of these cases, the balloon was placed adjacent to the aneurysm orifice, resulting in simultaneous **occlusion** of both the aneurysm and the parent artery, and in the other two cases, the balloon was positioned proximal to the aneurysm, resulting in temporary flow arrest. Three aneurysms were treated with either CAP or N-hexyl-**cyanoacrylate** after implantation of a balloon-expandable **tantalum** stent within the parent artery across the aneurysm orifice. Complete angiographic obliteration was achieved in all but one case. One aneurysm ruptured. Another partially **occluded** aneurysm reopened 10 weeks after treatment. In all cases treated without stents, distal migration of the polymer resulted in either stenosis or **occlusion** of the parent arteries. The combination of stent implantation and polymer injection resulted in permanent aneurysm **occlusion** without detectable polymer migration. An intravascular stent deployed within the parent artery across the aneurysm orifice acted as a safety net during endosaccular polymer injection by allowing blood to flow from the aneurysm cavity while preventing distal migration of liquid polymer.

CT Check Tags: Animal; Support, Non-U.S. Gov't
 Aneurysm: PA, pathology
 Aneurysm: RA, radiography
 *Aneurysm: TH, therapy
 Angiography, Digital Subtraction
 Balloon Dilatation
 Carotid Artery Diseases: PA, pathology
 Carotid Artery Diseases: RA, radiography
 *Carotid Artery Diseases: TH, therapy
 Cellulose: AD, administration & dosage
 *Cellulose: AA, analogs & derivatives
 Cellulose: TU, therapeutic use
 *Cyanoacrylates: AD, administration & dosage
 Cyanoacrylates: TU, therapeutic use
 Dogs
 Injections
 Polymers: AD, administration & dosage
 Polymers: TU, therapeutic use
 *Stents

RN 26877-39-4 (poly(hexyl-2-cyanoacrylate)); 9004-34-6 (Cellulose);
 9004-35-7 (acetylcellulose)

CN 0 (Cyanoacrylates); 0 (Polymers)

L15 ANSWER 12 OF 24 EMBASE COPYRIGHT 2002 ELSEVIER SCI. B.V.
 AN 95027311 EMBASE
 DN 1995027311
 TI Symptomatic pulmonary complications from liquid acrylate embolization of
 brain **arteriovenous malformations**.
 AU Pelz D.M.; Lownie S.P.; Fox A.J.; Hutton L.C.
 CS Radiology Department, University Hospital, Box 5339, 339 Windermere
 Rd, London, Ont. N6A 5A5, Canada
 SO American Journal of Neuroradiology, (1995) 16/1 (19-26).
 ISSN: 0195-6108 CODEN: AAJNDL
 CY United States
 DT Journal; Article
 FS 008 Neurology and Neurosurgery
 014 Radiology
 015 Chest Diseases, Thoracic Surgery and Tuberculosis
 037 Drug Literature Index
 038 Adverse Reactions Titles

LA English
 SL English
 TI Symptomatic pulmonary complications from liquid acrylate embolization of
 brain **arteriovenous malformations**.
 AB PURPOSE: To describe symptomatic pulmonary emboli from brain
arteriovenous malformation embolization with liquid
 acrylates and to analyze the reasons for these complications and describe
 preventive techniques. METHODS: The clinical records of 182 patients
 embolized with acrylate glue since 1978 for treatment of brain
AVMs were searched for evidence of symptomatic pulmonary
 complications. Originally isobutyl-2-cyanoacrylate and more
 recently n- butyl-2-cyanoacrylate were used in all patients.
Arteriovenous malformation morphology, amounts and
 techniques of glue injection, and clinical and radiologic investigations
 in the symptomatic patients were recorded. RESULTS: Three patients had
 pulmonary symptoms within 48 hours of glue injection. One patient with a
 left frontal **arteriovenous malformation** had
 embolization with an isobutyl-2-cyanoacrylate/pantopaque/acetic
 acid mixture; severe pleuritic chest pain developed 2 days later. One

patient with a left temporal and one with a left cerebellar **arteriovenous malformation** had embolization with n-butyl-2-cyanoacrylate/lipiodol mixtures; a cough, pleuritic chest pain, and bloody sputum developed in both within 24 hours. Two patients experienced a significant drop in PO2. No flow-arrest techniques were used for any of the injections in these three patients. All patients demonstrated significant changes on chest x-ray and CT chest examinations. All were treated conservatively and recovered spontaneously. **CONCLUSIONS:** Symptomatic pulmonary complications can occur after acrylate glue injection, particularly when delivery systems without flow arrest are used in high-flow vascular brain lesions. Techniques using acetic acid to delay polymerization time and 'sandwich' techniques in which glue is pushed with dextrose are also more susceptible to this complication.

CT Medical Descriptors:

*artificial embolism

***brain arteriovenous malformation: TH, therapy**

***brain arteriovenous malformation: DI, diagnosis**

*lung embolism: SI, side effect

adult

article

case report

clinical feature

computer assisted tomography

female

human

male

oxygen tension

thorax pain

thorax radiography

Drug Descriptors:

*acetic acid: CB, drug combination

*bucrilate: AE, adverse drug reaction

*bucrilate: CB, drug combination

***cyanoacrylate derivative: AE, adverse drug reaction**

*enbucrilate: AE, adverse drug reaction

*enbucrilate: CB, drug combination

*iodinated poppyseed oil: CB, drug combination

*iofendylate: CB, drug combination

glue

tantalum

RN (acetic acid) 127-08-2, 127-09-3, 64-19-7, 71-50-1; (bucrilate) 1069-55-2; (enbucrilate) 25154-80-7, 6606-65-1; (iodinated poppyseed oil) 8001-40-9, 8002-46-8, 8006-56-2, 8006-57-3; (iofendylate) 99-79-6; (**tantalum**) 7440-25-7

L15 ANSWER 13 OF 24 MEDLINE DUPLICATE 7

AN 92325777 MEDLINE

DN 92325777 PubMed ID: 1625008

TI **Arteriovenous malformations** of the brain: choosing embolic materials to enhance safety and ease of excision.

CM Comment in: J Neurosurg. 1993 Jul;79(1):153-5

AU Purdy P D; Batjer H H; Risser R C; Samson D

CS Department of Radiology, University of Texas Southwestern Medical School, Dallas.

SO JOURNAL OF NEUROSURGERY, (1992 Aug) 77 (2) 217-22.

Journal code: 0253357. ISSN: 0022-3085.

CY United States

DT Journal; Article; (JOURNAL ARTICLE)

LA English

FS Abridged Index Medicus Journals; Priority Journals
 EM 199208
 ED Entered STN: 19920821
 Last Updated on STN: 20000303
 Entered Medline: 19920813
 TI **Arteriovenous malformations** of the brain: choosing
 embolic materials to enhance safety and ease of excision.
 AB The authors report their experience with surgical resection of 108
 previously embolized **arteriovenous malformations** (**AVM**'s). Embolization was performed via only transfemoral
 catheterization in 70 lesions and via the surgical exposure of feeding
 vessels in 32. The remaining six patients were referred for resection
 following silicone sphere embolization elsewhere. Materials used included
 polyvinyl alcohol (PVA) foam, **platinum** microcoils, detachable
 silicone balloons, surgical silk, a mixture of 33% ethanol and
 microfibrillar collagen, and isobutyl **cyanoacrylate** (IBCA). It
 is believed that proximal arterial **occlusion** with balloons is an
 inferior choice for preresection embolization, because the technical
 difficulty of placement is high and the nidus of the **AVM** is
 unaffected. Vascular coagulation and section and **AVM** retraction
 are more difficult with IBCA; therefore, this is also considered an
 inferior choice. Among the materials studied, the combination of PVA for
 distal **occlusion** and microcoils for proximal **occlusion**
 appears to be the superior choice. Fewer complications (stroke or
 hemorrhage) are seen when intraarterial Amytal (amobarbital) testing is
 used to guide the embolization. Data regarding toxicity, oncogenicity, and
 vascular metabolism or recanalization associated with PVA, IBCA, and
 n-butyl **cyanoacrylate** are reviewed.
 CT Check Tags: Human
 Balloon Dilatation
 Bucrylate: PK, pharmacokinetics
 *Bucrylate: TU, therapeutic use
 *Embolization, Therapeutic
 Intracranial Arteriovenous Malformations: SU, surgery
 *Intracranial Arteriovenous Malformations: TH, therapy
 Polyvinyl Alcohol: PK, pharmacokinetics
 *Polyvinyl Alcohol: TU, therapeutic use
 L15 ANSWER 14 OF 24 MEDLINE DUPLICATE 8
 AN 85138293 MEDLINE
 DN 85138293 PubMed ID: 3974809
 TI Experimental carotid aneurysms: Part 2. Endovascular treatment with
cyanoacrylate.
 AU Kerber C W; Cromwell L D; Zanetti P H
 SO NEUROSURGERY, (1985 Jan) 16 (1) 13-7.
 Journal code: 7802914. ISSN: 0148-396X.
 CY United States
 DT Journal; Article; (JOURNAL ARTICLE)
 LA English
 FS Priority Journals
 EM 198503
 ED Entered STN: 19900320
 Last Updated on STN: 20000303
 Entered Medline: 19850326
 TI Experimental carotid aneurysms: Part 2. Endovascular treatment with
cyanoacrylate.
 AB Using our modification of the vein patch technique, we created 16
 aneurysms in the common carotid arteries of dogs. After a stabilizing and
 healing period, these aneurysms were treated using percutaneous catheter

techniques. Coaxial microcatheters were placed into the aneurysms, and a mixture of isobutyl 2-cyanoacrylate and tantalum dust was infused through the microcatheter using real time fluoroscopic control. Fifteen of the 16 aneurysms were successfully occluded; 1 was a failure because of total occlusion of the carotid artery. One human facial artery aneurysm was similarly treated. The ease and technical details of the treatment are discussed. Although the results are encouraging, we believe that it would be prudent to broaden the animal experimentation rather than to begin human use. Because no experimental aneurysm models are yet physiological, our results must be applied with caution to human intracranial aneurysms.

CT Check Tags: Animal; Case Report; Female; Human; Male

Adolescence

*Bucrylate

Carotid Artery Diseases: ET, etiology

*Carotid Artery Diseases: TH, therapy

Carotid Artery, External

*Cyanoacrylates

Dogs

*Embolization, Therapeutic: MT, methods

Intracranial Aneurysm: ET, etiology

*Intracranial Aneurysm: TH, therapy

Tantalum

RN 1069-55-2 (Bucrylate); 7440-25-7 (Tantalum)

CN 0 (Cyanoacrylates)

L15 ANSWER 15 OF 24 MEDLINE DUPLICATE 9

AN 82245789 MEDLINE

DN 82245789 PubMed ID: 7099369

TI The clinical application of intracranial artery cannulation technique (author's transl).

AU Negoro M; Berenstein A

SO NO SHINKEI GEKA. NEUROLOGICAL SURGERY, (1982 Mar) 10 (3) 271-7.

Journal code: 0377015. ISSN: 0301-2603.

CY Japan

DT Journal; Article; (JOURNAL ARTICLE)

LA Japanese

FS Priority Journals

EM 198209

ED Entered STN: 19900317

Last Updated on STN: 20000303

Entered Medline: 19820910

AB The introduction of the microscope into the neurosurgical operating theater brought the significant change on its operative results. However, even by means of the meticulous microsurgical techniques, certain intracranial lesion like deep-seated AVM cannot yet be successfully treated. Instead of the extravascular approach, intravascular treatment of these lesions has been evolved and become the great aid for the therapeutic purpose. In 1974 Serbinenko published his excellent work about his detachable balloon catheter technique. He succeeded in treating the intracranial lesions by the intravascular approach with the more exact manner than before. The balloon could make it possible to guide the small catheter into the distal branch of intracranial arteries. And also the balloon was detached and use as embolus. Until now, various balloon catheters become clinically available. Among them the catheter which Kerber devised is made of soft silicone and equipped with microballoon at the distal end. Although the balloon itself cannot be detached, it has a small hole at its top and can deliver the fluid through this opening (calibrated leak). The method for this catheter usage is as follows. Using

Seldinger technique, the non-tapered thin wall catheter has to be placed on the proximal side of the attempted artery as the introducing catheter. Through it, balloon catheter is cannulated coaxially and navigated more distally with the inflation or deflation of the balloon. Clinical application of this catheter include the superselective angiography, drug infusion and selective embolization. For the embolization, fluid embolus must be chosen. At this time **cyanoacrylate**, a potent tissue adhesive, is used as the embolus and injected with the mixture of Pantopaque and **tantalum** powder. Two cases of deep seated cerebral **AVM** were treated by selective embolization. **AVM** was completely **occluded** in one case, in the other case the embolization was interrupted because of worsening of neurologic deficits. In conclusion, the calibrated leak balloon catheter (Kerber) has wide range of clinically applicable potential and will become the great aid for the intravascular treatment.

CT Check Tags: Case Report; Human; Male

Adult

*Catheterization: MT, methods

Cerebral Angiography

*Cerebral Arteries: SU, surgery

*Embolization, Therapeutic: MT, methods

English Abstract

*Intracranial Arteriovenous Malformations: TH, therapy

Middle Age

L15 ANSWER 16 OF 24 EMBASE COPYRIGHT 2002 ELSEVIER SCI. B.V.

AN 82063297 EMBASE

DN 1982063297

TI [Superselective embolization in the genitourinary tract].
SUPERSELEKTIVER ARTERIELLER GEFASSVERSCHLUSS IM UROGENITALTRAKT.

AU Guenther R.

CS Inst. Klin. Strahlen., Johannes Gutenberg-Univ., 6500 Mainz 1, Germany

SO Aktuelle Urologie, (1982) 13/1 (1-4).

CODEN: AKURAJ

CY Germany

DT Journal

FS 028 Urology and Nephrology

018 Cardiovascular Diseases and Cardiovascular Surgery

014 Radiology

016 Cancer

037 Drug Literature Index

LA German

SL English

AB Superselective embolization of the kidney with pinpointed **occlusion** of small branches of the renal artery for treatment of a.-v. fistulae, bleeding due to angiomas or following biopsy is the method of choice with curative effect. In inoperable tumors of solitary kidneys as well as tumors of the bladder, prostate or uterus it is only a palliative measure generally with the aim of hemostasis. The tissue adhesive Butyl-2-**cyanoacrylate** mixed with lipiodol and **tantalum** powder has proved particularly suitable for superselective **occlusion** of small vessels.

CT Medical Descriptors:

*arteriovenous fistula

*artificial embolism

*hemangioma

*kidney artery

*kidney cancer

peripheral vascular system

therapy
kidney
Drug Descriptors:
 *cyanoacrylate derivative
 *gelfoam
 *enbucrilate
 *polyvinyl alcohol sponge
 *prolamin
 *silicone
 iodinated poppyseed oil
 tantalum
 ethibloc

RN (enbucrilate) 25154-80-7, 6606-65-1; (polyvinyl alcohol sponge)
63148-64-1; (prolamin) 117987-77-6; (silicone) 63148-53-8, 8043-93-4,
8055-24-1; (iodinated poppyseed oil) 8001-40-9, 8002-46-8, 8006-56-2,
8006-57-3; (**tantalum**) 7440-25-7; (ethibloc) 91196-33-7

L15 ANSWER 17 OF 24 EMBASE COPYRIGHT 2002 ELSEVIER SCI. B.V.

AN 82036447 EMBASE

DN 1982036447

TI Pathology of **arteriovenous malformations** embolized
with isobutyl-2-**cyanoacrylate** (bucrylate). Report of two cases.

AU Vinters H.V.; Debrun G.; Kaufmann J.C.E.; Drake C.G.

CS Dept. Clin. Neurol. Sci., Univ. Hosp., London, Ont. N6A 5A5, Canada

SO Journal of Neurosurgery, (1981) 55/5 (819-825).

CODEN: JONSAC

CY United States

DT Journal

FS 038 Adverse Reactions Titles

037 Drug Literature Index

008 Neurology and Neurosurgery

014 Radiology

005 General Pathology and Pathological Anatomy

018 Cardiovascular Diseases and Cardiovascular Surgery

LA English

TI Pathology of **arteriovenous malformations** embolized
with isobutyl-2-**cyanoacrylate** (bucrylate). Report of two cases.

AB There is controversy as to the possible toxic effects of isobutyl-2-**cyanoacrylate** (bucrylate) when this substance is used for purposes
of therapeutic embolization. Two cases are presented in which cerebral
arteriovenous malformations were resected, one 42 days
and the other a year after bucrylate embolization. In both, pathological
examination revealed a brisk intimal foreign-body giant-cell reaction
wherever bucrylate was present in a vessel, along with chronic
inflammation in the vessel walls and adjacent brain parenchyma. The
findings are discussed in the light of other observations on the
histotoxicity of bucrylate.

CT Medical Descriptors:

*adverse drug reaction

*artificial embolism

***brain arteriovenous malformation**

*embolism

*giant cell granuloma

*therapy

histology

central nervous system

peripheral vascular system

case report

autopsy

Drug Descriptors:

*bucrilate

*polyvinyl alcohol

***tantalum**

polyvinyl alcohol sponge

RN (bucrilate) 1069-55-2; (polyvinyl alcohol) 37380-95-3, 9002-89-5; (**tantalum**) 7440-25-7; (polyvinyl alcohol sponge) 63148-64-1

L15 ANSWER 18 OF 24 MEDLINE DUPLICATE 10

AN 81200154 MEDLINE

DN 81200154 PubMed ID: 6165036

TI Treatment of intracerebral **arteriovenous malformations** with isobutyl 2-**cyanoacrylate**: initial clinical experience.

AU Bank W O; Kerber C W; Cromwell L D

SO RADIOLOGY, (1981 Jun) 139 (3) 609-16.

Journal code: 0401260. ISSN: 0033-8419.

CY United States

DT Journal; Article; (JOURNAL ARTICLE)

LA English

FS Abridged Index Medicus Journals; Priority Journals

EM 198107

ED Entered STN: 19900316

Last Updated on STN: 19970203

Entered Medline: 19810723

TI Treatment of intracerebral **arteriovenous malformations** with isobutyl 2-**cyanoacrylate**: initial clinical experience.

AB From November 1976 to September 1979, 46 patients with intracranial **arteriovenous malformations** or fistulas participated in a clinical study using isobutyl 2-**cyanoacrylate** (IBCA), with **tantalum**, for palliative or preoperative **occlusion** of the blood supply to the abnormalities. Although failure to obtain satisfactory position of a functioning microcatheter precluded deposition of IBCA 10 times, a total of 51 of a possible 62 feeding vessels were **occluded** with the **tantalum**-impregnated glue. The technique, results, and complications are discussed in light of the clinical follow-up, which varied from 12 to 48 months.

CT Check Tags: Female; Human; Male

Adolescence

Adult

Aged

Arteriovenous Malformations: RA, radiography

***Arteriovenous Malformations: TH, therapy**

***Bucrylate: AD, administration & dosage**

Child

***Cyanoacrylates: AD, administration & dosage**

***Embolization, Therapeutic: MT, methods**

Follow-Up Studies

Middle Age

Palliative Care

Postoperative Complications

Tantalum

Tomography, X-Ray Computed

RN 1069-55-2 (Bucrylate); 7440-25-7 (**Tantalum**)

CN 0 (**Cyanoacrylates**)

L15 ANSWER 19 OF 24 EMBASE COPYRIGHT 2002 ELSEVIER SCI. B.V.

AN 81070275 EMBASE

DN 1981070275

TI [Transrenal ureteric embolisation. Experimental and clinical results].

TRANSRENALE URETEREMBOLISATION. EXPERIMENTELLE UND KLINISCHE ERGEBNISSE.

AU Guenther R.; Klose K.; Bohl J.; Marberger M.
 CS Inst. Klin. Strahlenk., Univ. Mainz, Germany
 SO Fortschritte auf den Gebiete der Rontgenstrahlen und der Nuklearmedizin,
 (1980) 133/5 (471-476).
 CODEN: FGRNAJ

CY Germany
 DT Journal
 FS 014 Radiology
 028 Urology and Nephrology
 037 Drug Literature Index

LA German
 SL English

AB Transrenal ureteric embolisation with the tissue adhesive butyl
 2-cyano-acrylate mixed with lipiodol and **tantalum** powder
 produces rapid and effective **occlusion** of the ureter. This can
 be combined with the use of a Gianturco spiral. In an experimental series
 of eight dogs observed for 54 days, the adhesive disappeared partly in
 three animals and totally in a further three. Fibrotic ureteric stenosis
 was observed in three animals. Clinical results were rather better because
 of the presence of external urinary drainage. Twenty-one ureters in
 eighteen patients were **occluded** either unilaterally or
 bilaterally, and sometimes combined with percutaneous contralateral renal
 embolisation. After one to four months, the ureter was still
occluded in three out of six cases; after five to 17 months, three
 out of five cases were no longer totally **occluded**, but in two
 cases they were still blocked. The procedure is new and suitable for the
 treatment of otherwise untreatable conditions, such as extensive urinary
 fistulae, bladder tenesmus, and haematuria due to extensive tumours of the
 minor pelvis.

CT Medical Descriptors:
 *artificial embolism
 *ureter
 bladder carcinoma
 uterine cervix carcinoma
 urinary tract
 bladder
 animal experiment
 major clinical study
 therapy
 dog
 kidney
 Drug Descriptors:
 *cyanoacrylate derivative
 *iodinated poppyseed oil
 *tantalum

RN (iodinated poppyseed oil) 8001-40-9, 8002-46-8, 8006-56-2, 8006-57-3; (
tantalum) 7440-25-7

L15 ANSWER 20 OF 24 EMBASE COPYRIGHT 2002 ELSEVIER SCI. B.V.
 AN 79176753 EMBASE
 DN 1979176753
 TI Modification of **cyanoacrylate** for therapeutic embolization:
 Preliminary experience.
 AU Cromwell L.D.; Kerber C.W.
 CS Dept. Radiol., Univ. Washington Sch. Med., Seattle, Wash. 98195, United
 States
 SO American Journal of Roentgenology, (1979) 132/5 (799-801).
 CODEN: AJROAM

CY United States
 DT Journal
 FS 037 Drug Literature Index
 014 Radiology
 018 Cardiovascular Diseases and Cardiovascular Surgery
 LA English
 TI Modification of **cyanoacrylate** for therapeutic embolization:
 Preliminary experience.
 AB **Cyanoacrylate** is a satisfactory material for therapeutic
 embolization, but it has the disadvantages of not being radiopaque and
 polymerizing within 1 sec after contact with ionic materials. Its behavior
 was modified with varying concentrations of iophendylate and we were able
 to satisfactorily control its polymerization from 1 to 30 sec. This
 control should allow penetration of **arteriovenous**
malformations, which is necessary if cure is to result. The
 iophendylate adds radiopacity and seems to enhance the suspension of
tantalum, another opacifying agent. Preliminary experience in dogs
 is encouraging, but too few humans have been treated with this method to
 recommend it as more than an experimental procedure at this time.
 CT Medical Descriptors:
 *artificial embolism
 dog
 methodology
 therapy
 animal experiment
 peripheral vascular system
 Drug Descriptors:
 ***cyanoacrylate**
 *iophendylate
 RN (**cyanoacrylate**) 15802-18-3; (iophendylate) 99-79-6
 L15 ANSWER 21 OF 24 MEDLINE DUPLICATE 11
 AN 79247739 MEDLINE
 DN 79247739 PubMed ID: 472240
 TI Catheter and material selection for transarterial embolization: technical
 considerations. II. Materials.
 AU Berenstein A; Kricheff I I
 SO RADIOLOGY, (1979 Sep) 132 (3) 631-9.
 Journal code: 0401260. ISSN: 0033-8419.
 CY United States
 DT Journal; Article; (JOURNAL ARTICLE)
 LA English
 FS Abridged Index Medicus Journals; Priority Journals
 EM 197910
 ED Entered STN: 19900315
 Last Updated on STN: 19970203
 Entered Medline: 19791017
 AB In this second part of the report, the authors discuss the advantages and
 disadvantages of several embolization agents. These include Gelfoam,
 silicone spheres, polyvinyl alcohol foam (PVA), isobutyl-2-
cyanoacrylate (IBCA), silicone fluid mixtures, and
tantalum powder. The techniques employed and conditions under
 which these materials should be used are discussed.
 CT Check Tags: Animal; Comparative Study; Female; Human; Male
 Adult
 Aged
 Angiography
 Arteriovenous Malformations: RA, radiography
 Arteriovenous Malformations: TH, therapy

Bucrylate
Carcinoma, Squamous Cell: BS, blood supply
Carcinoma, Squamous Cell: RA, radiography
Carcinoma, Squamous Cell: TH, therapy
Catheterization: IS, instrumentation
Child, Preschool
Embolization, Therapeutic: IS, instrumentation
*Embolization, Therapeutic: MT, methods
Facial Neoplasms: RA, radiography
Facial Neoplasms: TH, therapy
Gelatin Sponge, Absorbable
Meningioma: BS, blood supply
Meningioma: RA, radiography
Meningioma: TH, therapy
Microspheres
Middle Age
Polyvinyl Alcohol
Silicones

Tantalum

RN 1069-55-2 (Bucrylate); 7440-25-7 (Tantalum); 9002-89-5
(Polyvinyl Alcohol)

L15 ANSWER 22 OF 24 EMBASE COPYRIGHT 2002 ELSEVIER SCI. B.V.

AN 79150782 EMBASE

DN 1979150782

TI **Cyanoacrylate occlusion** of carotid-cavernous fistula
with preservation of carotid artery flow.

AU Kerber C.W.; Bank W.O.; Cromwell L.D.

CS Dept. Radiol., Univ. Pittsburgh Sch. Med., Pittsburgh, Pa. 15261, United
States

SO Neurosurgery, (1979) 4/3 (210-215).

CODEN: NRSRDY

CY United States

DT Journal

FS 037 Drug Literature Index

008 Neurology and Neurosurgery

018 Cardiovascular Diseases and Cardiovascular Surgery

LA English

TI **Cyanoacrylate occlusion** of carotid-cavernous fistula
with preservation of carotid artery flow.

AB We report a new treatment for carotid-cavernous fistula. Using a flow-guided, balloon-tipped microcatheter we catheterize the fistula itself, verify balloon entry into the fistula with fluoroscopy and x-ray films, and then infuse the tissue adhesive isobutyl-2-**cyanoacrylate** with careful fluoroscopic control. Three patients have had their fistulas **occluded**, with preservation of flow through the internal carotid artery. This balloon microcatheter allows the radiologist to perform a reversible test **occlusion**. All three patients had neurological changes during or after the procedure, and in one we inadvertently **occluded** several distal middle cerebral artery branches without permanent neurological deficit. No patient became blind or developed 3rd, 4th, or 6th nerve palsy from the treatment. This technique seems to have promise as another method for the obliteration of carotid-cavernous fistula.

CT Medical Descriptors:

*artificial embolism

*carotid cavernous fistula

*carotid artery

*carotid artery fistula

*carotid artery flow
*cavernous sinus carotid artery fistula
balloon catheter
drug therapy
peripheral vascular system
major clinical study
therapy
Drug Descriptors:
*bucrilate
 ***cyanoacrylate**
 ***tantalum**

diazepam

RN (bucrilate) 1069-55-2; (**cyanoacrylate**) 15802-18-3; (**tantalum**) 7440-25-7; (diazepam) 439-14-5

L15 ANSWER 23 OF 24 MEDLINE

AN 79125916 MEDLINE

DN 79125916 PubMed ID: 570455

TI Transcatheter embolization of the kidney with butyl-2-**cyanoacrylate**: experimental and clinical results.

AU Gunther R; Schubert U; Bohl J; Georgi M; Marberger M

SO CARDIOVASCULAR RADIOLOGY, (1978 Apr 25) 1 (2) 101-8.

Journal code: 7807044. ISSN: 0342-7196.

CY United States

DT Journal; Article; (JOURNAL ARTICLE)

LA English

FS Priority Journals

EM 197905

ED Entered STN: 19900315

Last Updated on STN: 19970203

Entered Medline: 19790524

TI Transcatheter embolization of the kidney with butyl-2-**cyanoacrylate**: experimental and clinical results.

AB The technique and efficacy of therapeutic catheter embolization of the kidney with butyl-2-**cyanoacrylate** (Histoacryl) were studied in 80 rabbits (including control groups) and in 10 dogs. A mixture of butyl-2-**cyanoacrylate**, 50% glucose, and **tantalum** powder was used for the embolization. Complete and permanent vascular **occlusion** was found in nearly all cases. The main complication observed was a reflux of embolizing material into the lumbar arteries, which occurred in seven rabbits. Clinically therapeutic embolization was performed in six patients with hypernephroma. The indication for embolization in these patients, as well as in two others with iatrogenic lesions, was pronounced hematuria. Cessation of bleeding was achieved in all cases. For embolization the coaxial catheter technique is recommended; in special cases with extensive arteriovenous shunts, adjunctive balloon **occlusion** would be advisable.

CT Check Tags: Animal; Human

Adenocarcinoma: CO, complications

Catheterization: MT, methods

***Cyanoacrylates: TU, therapeutic use**

Dogs

Embolization, Therapeutic: AE, adverse effects

***Embolization, Therapeutic: MT, methods**

***Enbucrilate: TU, therapeutic use**

Hematuria: ET, etiology

***Hematuria: TH, therapy**

Kidney Neoplasms: CO, complications

Lumbosacral Region: BS, blood supply

Rabbits

*Renal Artery

Renal Artery: RA, radiography

Tantalum: TU, therapeutic use

RN 6606-65-1 (Enbucrilate); **7440-25-7 (Tantalum)**

CN 0 (**Cyanoacrylates**)

L15 ANSWER 24 OF 24 EMBASE COPYRIGHT 2002 ELSEVIER SCI. B.V.

AN 78231800 EMBASE

DN 1978231800

TI [The technique of therapeutic catheter embolization of the kidney with Histoacryl (Butyl 2 **cyanoacrylate**)].

TECHNIK DER THERAPEUTISCHEN KATHETEREMBOLISIERUNG DER NIERE MIT HISTOACRYL (BUTYL 2 CYANOACRYLAT).

AU Guenther R.; Schubert U.; Georgi M.; Marberger M.

CS Inst. Klin. Strahlenkunde, Univ. Mainz, Germany

SO Aktuelle Urologie, (1977) 8/6 (229-303).

CODEN: AKURAJ

CY Germany

DT Journal

FS 037 Drug Literature Index

028 Urology and Nephrology

014 Radiology

016 Cancer

LA German

SL English

TI [The technique of therapeutic catheter embolization of the kidney with Histoacryl (Butyl 2 **cyanoacrylate**)].

TECHNIK DER THERAPEUTISCHEN KATHETEREMBOLISIERUNG DER NIERE MIT HISTOACRYL (BUTYL 2 CYANOACRYLAT).

AB On account of experimental research of 50 rabbits and 10 dogs and based on empirical data of 7 patients the authors are able to give exact details of a practical and successful application of Histoacryl (butyl-2-cyanoacrylat) for therapeutic catheter embolization of the kidney.

Histoacryl has to be diluted with glucose solution in an adequate ratio of components and made radiopaque by added **tantalum** powder. For the embolization it is advisable to use the coaxial catheter technique and the adjuvant balloon **occlusion** for special cases.

CT Medical Descriptors:

*artificial embolism

*catheter

*dog

*kidney cancer

*rabbit

methodology


theoretical study

Drug Descriptors:

***cyanoacrylate derivative**

enbucrilate

RN 178328-40-0 REGISTRY
CN 2-Propenoic acid, 3-cyano-, 6-cyanoheptyl ester (9CI) (CA INDEX
NAME)
OTHER NAMES:
CN 6-Cyanoheptyl cyanoacrylate
FS 3D CONCORD
MF C11 H14 N2 O2
SR CA
LC STN Files: CA, CAPLUS

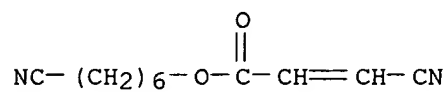

PROPERTY DATA AVAILABLE IN THE 'PROP' FORMAT